




Zitzi
SEABASS PRO



Crashtested and approved!

Seabass PRO is crashtested and approved according to ISO 7176-19 for 135 kg and Hi-tested according to ISO 10542.

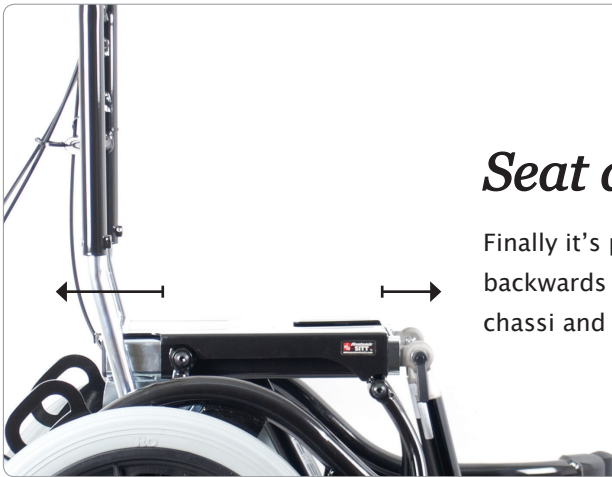

Seat depth adjustment

Finally it's possible to adjust the seat depth both forwards and backwards so that you can easily find the perfect balance between this chassi and any number of seating systems and/or custom made seats.



Tilt in space

Unique for Seabass PRO is its "tilt in space" feature, allowing the seat to be tilted without compromising the balance and stability between the seat and chassi. This way the center of gravity remains unchanged so that the seat system is always perfectly balanced regardless of the degree of tilt.



Wheels

In order to fulfill our customers' wishes, rear wheels are available in 5 different sizes. Choose between 12", 16", 20", 22", and 24" wheels. All of the rear wheels are equipped with quick-release function and are available with either air-filled or solid tires. Choose the optional e-drive and roll yourself or run the engine.

The picture is showing
A Seabass Pro size 1 with 16" wheels

Adjustable driving bow
Unique design for comfortable handling



e-drive
E-drive package
A flexible way to transfer

Flat seat plate
Pre-made tracks for an easier mount of adapter systems

Different wheels
To meet many needs

Abduction frame
To avoid problems between the footrest and casters

Quick-release axle
Easy to release and re-mount

Transportation kit
Approves transportation according to ISO 7176-19

Seabass Pro Quickies

Seabass Pro with its center of gravity on the right place. Versatile wheel location for a huge variety and wheel sizes. Whatever seating solution chosen this guarantees perfect balance for every fitting.

Sizes..... 3 st
Wheels 12-24"
Load weight 135 kg
Tilt 35°

Seat height.....41/46 cm
Color.....Black
Weight20 kg
Crash-testedApproved