



easicare

making things better

Thunder



Compact Mid-Wheel Drive Power Chair

The Easicare Thunder is a class leading compact mid-wheel drive power chair, with an inherent level of adjustability and comfort. The Thunder is built of high quality components such as the Dynamics Electronic Control System to ensure long term reliability and yet, at the same time, providing exceptional value.

Superb Indoor Manoeuvrability

The Easicare Thunder is truly nimble with a compact size and mid-wheel drive system, combining for a tight turning radius of only 620mm / 25 inches, making indoor use safer and easier.

Superb Outdoor Performance

With its twin motor drive system, a top speed of 6.7Kmh / 4.2mph and a travel range of 20Km / 12.4 miles, the Thunder provides superb outdoor performance. The mid-wheel drive system utilises 10" drive wheels, complimented by 6" castor wheels front and rear, to give controlled safety and stability over uneven surfaces.

Mechanical Adjustability

Inherent modular adjustability is included in the design specification resulting in easy and quick to adjust settings to accommodate different user needs.



Footrest in the upright position




Armrests detach for ease of transfer




Backrest folds for ease of storage



Positioning belt as standard

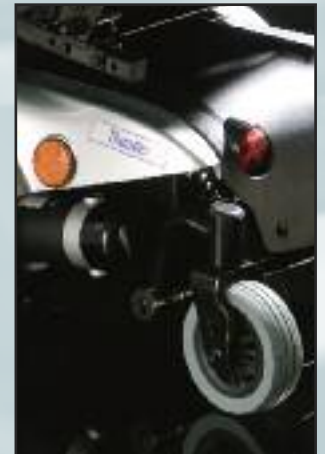

Light reflectors, castor and suspension

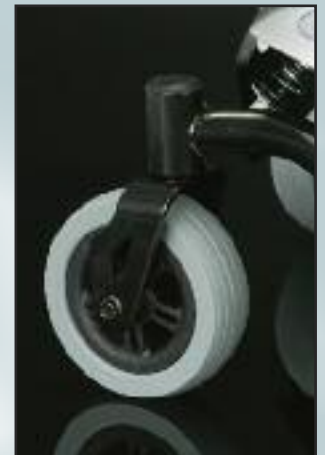

Front castor



Controller



Height adjustable armrest



Adjustable armrest



Full coil-over suspension system



Full swivel seat system



Adjustable backrest angle from 100° - 120°



Seat swivelled



Easy to reach charging socket



AVAILABLE OPTIONS:

- Transport headrest
- Tray
- Attendant control
- Swing away controller
- 40Ah battery upgrade

FEATURES:

- Maximum Speed: 6.7kmph / 4.2mph
- Travel Range: 20km / 12 miles (36Ah);
25km / 15.5 miles (40Ah);
- Capacity: 136Kg / 300lbs
- Full Coil - Over Suspension
- Adjustable Seat Width / Depth
- Adjustable Backrest Angle
- Height Adjustable Arm Rest
- Left / Right Adjustable Controller
- Length Adjustable Controller

SPECIFICATIONS:

- Overall Length: 995mm / 39"
- Overall Width: 610mm / 24"
- Overall Height: 930mm / 36.6"
- Wheels - Front: 155mm / 6"
- Wheels - Middle: 260mm / 10"
- Wheels - Rear: 155mm / 6"
- Weight with Batteries: 80Kg / 176lbs
- Ground Clearance: 60mm / 2.4"
- Grade Climbable: 6 degrees
- Kerb Climbing: 40mm / 1.6"
- Turning Radius: 620mm / 24.4"
- Brake: Electro-Mechanical
- Seat Width: 405 - 510mm / 16" - 20"
- Seat Depth: 405 - 460mm / 16" - 18"
- Seat Height: 430 - 480mm / 17" - 19"
- Back Seat Angle: 100 degree - 120 degrees
- Back Seat Height: 460mm / 18"
- Drive Train: 2-Motor mid-wheel drive
- Battery Weight: 24.2Kg / 53lbs
- Motor Size: 420W 4900rpm
- Battery: (2) 12V. 36Ah / (2) 12V. 40Ah
- Charger: 5A Off Board
- Electronics: Dynamics Electronic System
- Seat Type: Padded, breathable ballistic nylon back rest and positioning belt



www.mmsmedical.ie

The Easicare Thunder Power Chair is available from

MMS Medical.
51 East Gate Drive,
Little Island
Cork.

T: 021 4618000
F: 021 4618099
E: info@mmsmedical.ie

MMS Medical
Unit 81,
Western Parkway Business Park,
Lower Ballymount Road,
Dublin 2

T: 01 4500582
E: sales@mmsmedical.ie

For further information on the full specification of options and accessories of the Easicare Thunder Power Chair, please refer to the order form. All information is subject to change without notification. Please consult your dealer with any queries you may have. E&OE.