



MMS
MEDICAL

Twin rail curved stairlifts

Handicare 2000 | 2000 Simplicity



handicare

The benefits of a stairlift from Handicare

If you have difficulty getting up and down the stairs you are not alone. Every year thousands of people look for a solution to the problems posed by the difficulty they have using the staircase in their home.

The decision to install a stairlift is sometimes a difficult one to take but when compared with the alternative of moving house, or living on one level of their house, many people have found it to be their preferred option.

If you are thinking about installing a stairlift, and your stairs turn corners or are spiral, the Handicare 2000 series offers you a safe and reliable solution.


Handicare's stairlifts have helped thousands of people overcome the problems posed by going up and down their stairs at home.

The Handicare 2000 series stairlift suits most staircases with bends or corners and offers you all the basic functions that you would expect plus the option to upgrade to a folding footplate and/or powered seat swivel.



A tailor-made solution

Perfectly finished, 100% made to measure... and yet rapidly installed at your home. We will, after an informal visit, explain all the features of our stairlifts to you in detail.

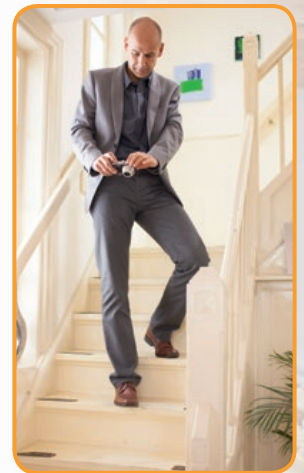

Each Handicare stairlift is designed to suit your individual staircase. We use advanced PhotoSurvey technology to make accurate measurements and we have all the options available to ensure the best possible fit for your staircase.

Your Handicare 2000 stairlift can be fitted with two seats. Your Handicare supplier will be able to advise you which seat is most suited to your needs.



Smart Seat



Simplicity Seat



Seat features



Remote control

The Handicare 2000 series comes with two remote controls which means, once folded, you can easily send it to the upstairs charge point freeing up space in your hallway.



Safety as standard

Your Handicare 2000 series stairlift will come with a range of safety sensors that are designed to stop your stairlift if it meets any obstruction as you travel up the stairs.



Foldable

The compact unit folds up leaving the staircase free for other users and the bespoke track is designed to fit your individual staircase and maximise clear space on the stairway. It can be supplied with a powered hinged track if required.

Battery

Like all Handicare stairlifts the 2000 series have a battery backup which means it will continue to work in the event of a power out so you do not need to worry about being stranded up- or downstairs.





Easy to use

The operating toggle is easy to use, simply hold it towards the direction you wish to travel, up or down, and your stairlift will move smoothly and gently on its way.



Seat belt

Your Handicare 2000 series have an easy release retractable seat belt for added peace of mind. If you prefer you can have an easy fasten Velcro seat belt.



Diagnostic display

The Handicare 2000 series stairlifts have a diagnostic display so that you can monitor the status of your stairlift, e.g. a safety edge has been activated.



Comfort features

Powered hinge

Some houses will have a doorway close to the bottom of the stairs and the track may intrude across this. In these circumstances an automatic hinged track should be fitted. This means that the track can be folded away from the doorway giving clear access. When operating the stairlift from the toggle, the hinge will lift/lower automatically.



Powered swivel

When you arrive at the top of the stairs you will need to swivel your seat round so that you exit onto the landing safely. If you choose the powered swivel option the swivel is operated by continuing to hold the toggle when you arrive at the top of the staircase. The seat will move automatically.

Weight capacity

The standard maximum weight limit on the Handicare 2000 series is 115kg this weight limit can be increased to 137kg by having the heavy duty kit fitted on your stairlift.

Powered footplate

The footplate and seat can be folded up to increase available space leaving the staircase as clear as possible for other users. If you choose the powered footplate option it will remove the need for you to bend down to fold or unfold the footplate.



Stairlift colours

Fitted to the stairs, not the wall, it is unlikely you will need to redecorate following the installation of your Handicare 2000 series stairlift. With a choice of six upholstery colours on the Smart seat and three track colours, your stairlift will compliment most colour schemes.

Perch

If you have restricted movement in the knee or hip joint you may find sitting painful. In these situations a Perch seat may be the solution.

A Perch seat means that the user is almost in a standing position when they are using the lift. There is a shortened seat and retractable seat belt which provide a stable support area for users to 'perch' on. Upholstery is available in a choice of six colours.




Handicare 2000 series: at a glance

Seat options	2000	2000 Simplicity	
Choice of Upholstery	✓	-	A selection of six colours means you can choose upholstery to match your décor.
Powered Swivel	○	-	This makes turning the stairlift seat to get in and out on the landing as easy as possible.
Retractable Seat Belt	✓	✓	This retractable belt can be clipped together on one side or centrally.
Velcro Seat Belt	○	○	This black belt has velcro on either end to hold it firmly in place.
Perch Seat	○	-	Ideal for users who have difficulty bending at the hip or knee joint.
Footplate options	2000	2000 Simplicity	
Manual	✓	✓	The footplate is folded and unfolded manually.
Powered from Seat	○	○	Folds the footplate away when the seat is folded up.
Powered from Arm	○	-	The footplate can be lifted and lowered by a switch under the arm whilst the user is seated.
Track Features	2000	2000 Simplicity	
Travels around corners and bends	✓	✓	For a wide range of staircases, from ones that turn a corner to multi-flight or spiral.
Multiple charge points	✓	✓	Use your remote control to park on charge points.
Track available in any colour	○	○	Customise your curved stairlift track to match your décor.
Zero intrusion into landing	○	○	For use when a doorway close to the top of the stairs would be obstructed by the track.
Powered Hinge	○	○	Automatically moves the track while the stairlift moves leaving doorways unobstructed.
Weight capacity	2000	2000 Simplicity	
Standard Maximum Weight	115 kg	115 kg	
Max weight with Heavy Duty Kit	137 kg	137 kg	
Standards compliance	2000	2000 Simplicity	
EN 81-40: Stairlifts	✓	✓	


✓ Available as standard ○ Optional - Not available

Upholstery Colour Choices Smart Seat

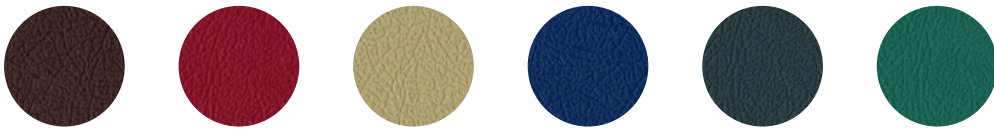

COCOA RUBY SAND SAPPHIRE SLATE JADE


Upholstery Colour Simplicity Seat



SAND

Please be aware that actual colours may vary slightly from those shown.

Track Colour Choices 2000 Series



LIGHT GREY BROWN GUN METAL

Handicare also offers special orders so the track can be finished in any colour to suit your décor.

Making everyday life easier

Experience

Convenience serves people... and lifting equipment of Handicare has been doing just that for more than 125 years.

With technical perfection, intelligent solutions and safe applications, Handicare's innovations have provided not only greater comfort, but also more attractive stairlifts... in your style, your taste.

Top quality in all respects make the Handicare stairlift a valuable asset for independence around the home. Many satisfied users experience that every day.

Quality and safety

All our new products are delivered with 2 years warranty.

A Handicare stairlift stands for high quality: quality in terms of safety, innovation, comfort and design. Furthermore, all our products are tested and approved by worldwide independent inspection services. This gives you the assurance that all our products meet the highest standards in terms of quality and safety.



MMS Medical,
51 East Gate Drive,
Little Island,
Cork,

T: 021 4618000
F: 021 4618099
E: info@mmsmedical.ie

MMS Medical,
Unit 81, Western Parkway Business Park,
Lower Ballymount Road,
Dublin 2

T: 01 4500582
F: 01 4565889
E: sales@mmsmedical.ie

www.mmsmedical.ie
www.stair-lifts.ie



handicare

Handicare Stairlifts B.V.
Newtonstraat 35
1704 SB Heerhugowaard
The Netherlands

T +31 (0)72 576 88 88
F +31 (0)72 574 34 35
stairlifts@handicare.com
www.handicare.com

Please be aware that actual colours may vary slightly from those shown. All information is provided in good faith but does not constitute a guarantee. Handicare reserves the right to alter designs and specifications without prior notice.