



Terry *Harmony*[®]

Home Lift



easy INDEPENDENT LIVING

Regain the freedom of your home with the Harmony Home Lift.

The lift will provide you with safe and easy access to your upstairs, meaning you can regain the whole living space without needing to move to a bungalow or ground floor flat.

Whether travelling seated in a wheelchair or on the optional fold-down seat, the Harmony Home Lift will allow you to move between floors safely, smoothly and quietly. The lift features an exceptionally wide range of car sizes to choose from, so you are sure to find a model that fits your needs. The design of the lift maximises the space in the car so you can travel in comfort.

The Harmony Lift can be fitted almost anywhere in your home. It requires only minimal headroom on the first floor and it can be installed against non-load bearing walls.

WE WILL TAKE CARE OF EVERYTHING

Purchasing and having your Harmony Lift installed is far less complicated than you might think.

We will survey your home and provide you with a written quotation. After your decision to proceed, we will first carry out all the necessary work. The installation itself takes less than two days.

PRODUCT *highlights*

- ✓ **No lift shaft required** - minimises the space taken up in your home
- ✓ **A market-leading selection of car sizes** to suit every need
- ✓ **Smart car design** minimises the lift's footprint
- ✓ **One-touch controls** – no need for constant pressure, like many low speed lifts
- ✓ **Low floor** ensures easy wheelchair access & egress
- ✓ **Quick installation** – fitted in less than two days
- ✓ Designed & manufactured **in the UK**

peace OF MIND

The safety of you and your family is our priority. This is why the Harmony Home Lift is equipped with all the essential safety features as standard.

Pressure sensitive safety surfaces stop the lift if it gets obstructed from above or below. The lift also features a battery back-up system, allowing you to return to the ground floor safely in the event of a power failure.

Fire safety is an important consideration. The Harmony Home Lift maintains fire and smoke protection between the floors, whether it is parked upstairs or downstairs. Smoke and heat alarms ensure the lift performs correctly in case of activation.

THE TERRY LIFTS *guarantee*

We guarantee the quality and safety of your home lift.

- ✓ Pressure sensitive surfaces and safe edges
- ✓ Emergency lowering in case of power cuts
- ✓ 30-minute fire and smoke integrity maintained between floors
- ✓ Smoke/heat alarms provided on both floors helping to ensure your safety
- ✓ Complies with the Machinery Directive and we are certificated by an Independent Notified Body for a lift travel above 3 metres (a legal requirement)
- ✓ Compliant with British Standard BS 5900:2012
- ✓ Fully compliant as standard with the Chartered Institute of Building Services Engineers (CIBSE) Guide D recommendations
- ✓ 2-year warranty including two service visits



comfortable AND EASY TO USE

Using your Harmony Home Lift is easy as 1, 2, 3.



Easy-to-use wireless controls come as standard in the car and at both floors.



A simple push of the button sends the lift up. The controls can be positioned to suit your needs.



The aperture ceiling panel is automatically lifted or replaced en-route.



The lift moves to the other floor and vacates the space it was in.



Low car construction ensures easy wheelchair access and egress.



The lift can be sent back to the floor below, leaving a feeling of space in the room.

your HARMONY LIFT

Standard features

Safe edges also stop the lift from moving if obstructed

Toughened safety glass windows, which do not scratch

Modern LED lighting illuminates the car and the room

Pressure-sensitive safety surface stops the lift when obstructed. The infill provides a minimum of 30-minute fire protection and seals the room from smoke

Lifting mechanism is concealed in the guide, resulting in a neat and tidy appearance

Adjustable, wireless car controls can be positioned on the left or right

Powered door is available on the wheelchair models

Pressure-sensitive safety surface, doubling as ceiling closure and fire/smoke seal


Optional extras

- ✓ A range of trim colours
- ✓ Folding seat (can be used as a perch seat)
- ✓ Telephone in the lift car
- ✓ Hand held remote control
- ✓ Grab handles
- ✓ Joystick control
- ✓ Left or right door handing
- ✓ Car padding on three sides and door
- ✓ Mirror
- ✓ Up to 5-year warranty

Harmony S

Car details

All dimensions in mm

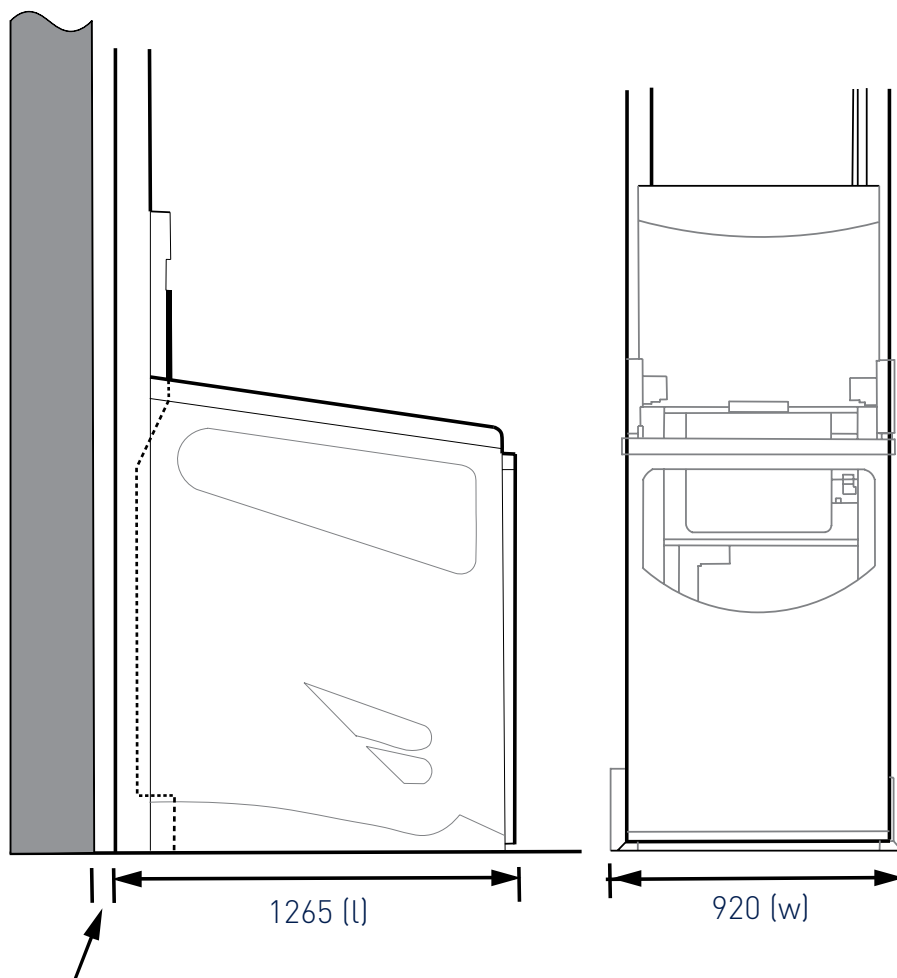

50 minimum
recommended

Travel: From 2060 to 3600

Lift Model

| Lift Model | Internal Width | Internal Length | External Width (w) | External Length (l) | Capacity |
|----------------------------|----------------|-----------------|--------------------|---------------------|----------|
| Harmony S (Std wheelchair) | 755 | 1175 | 920 | 1265 | 280 kg |
| Harmony W (Wider) | 890 | 1175 | 1055 | 1265 | 280 kg |
| Harmony L (Longer) | 755 | 1375 | 920 | 1465 | 250 kg |
| Harmony LW (Longer Wider) | 890 | 1375 | 1055 | 1465 | 250 kg |



For further information or to book your free home survey, call MMS Medical today on **021 4618000** or visit www.mmsmedical.ie to view our full range of home lifts



Lifestyle Home Lift



Vivendi Home Lift



Terry *Harmony*[®]

MMS Medical,
51 Eastgate Drive,
Little Island,
Cork
T: 021 4618000
E: info@mmsmedical.ie

