

QUICKIE®

A True Paediatric Chair that Grows with the Child



The Quickie RX Kidz - the First Choice Chair for Paediatric Users

Born from the Quickie RXS, the RX Kidz has been designed with the same great features and functionality, but is specifically for paediatric users.

Outrigged for Stability


The Quickie RX Kidz is available with swing-away or outrigged front frame. The latter increases forward stability, allowing the child to gain confidence without the risk of falls when leaning forward or transferring from or to the chair.

Room for Growth

The Quickie RX Kidz can grow with the child. It features infinitely adjustable legrests and footplates and an optional 2.5cm or 5cm seat depth growth kit.




Technical Specifications

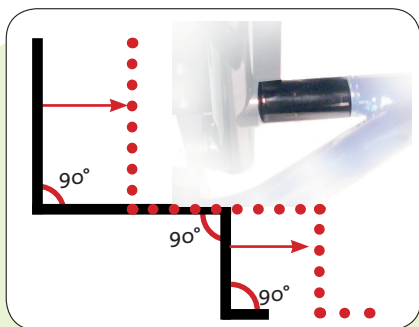

Optional 2.5cm or 5cm seat depth growth kit. 90° x 90° x 90° seating position always maintained through seat and frame growth kit and adjustable footplates



Adjustable legrests provide 6" of depth and 9" of height adjustment. No dead spots - infinite adjustment.

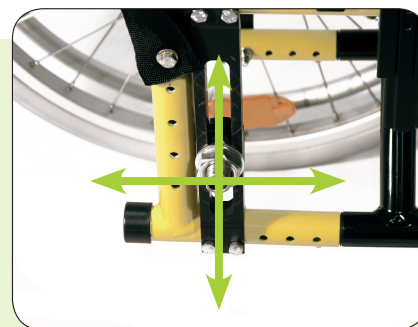

Centre of gravity and seat height adjustment via vertical axle plate (option)



Optional Squiggle upholstery and transit approved removable domed headrest.



Outriggered frame option for increased stability.



Crash tested.

Technical Specification

Seat Width:	25-38cm (10-15")	Front Seat Height:	44 - 51cm (17" - 20")
Seat Depth:	25-38cm (10-15")	Rear Seat Height:	36 - 49cm (14" - 19")
Back Height:	30.5-43cm (12-17")	Overall Width:	SW +22.5cm (9")
Back Angle:	8° half folding	Chair Weight:	13.4 kg
Frame Angles:	70° swing-away or 90° outriggered	Max User Weight:	75 kg (11.8 stone)
Colours:	Black, yellow, purple, red, blue, green		

For further information on the full specification, options and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.