

The background of the advertisement features four red dice on a white surface. One die is in sharp focus in the foreground, showing the 1, 2, and 3 faces. Three other dice are in the background, slightly out of focus, showing various faces. The overall lighting is bright and even.

BREEZY[®] BasiX[™] RubiX[™]

Play the
Combinations


Numerous options,
multiple uses,
one family!







Play the Combinations!

- ONE common system for adult wheelchair applications
 - BasiX for part time and passive use
 - RubiX for full time, hemiplegic, bariatric, semi active and passive use
- Easily adjustable seat depth, height and angle
- Interchangeable options and accessories for easy conversion
- Durable 'all aluminium' design for long life
- Fully transit tested
- Easy to order

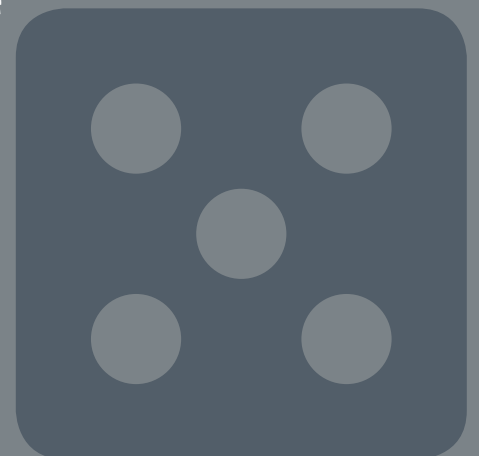

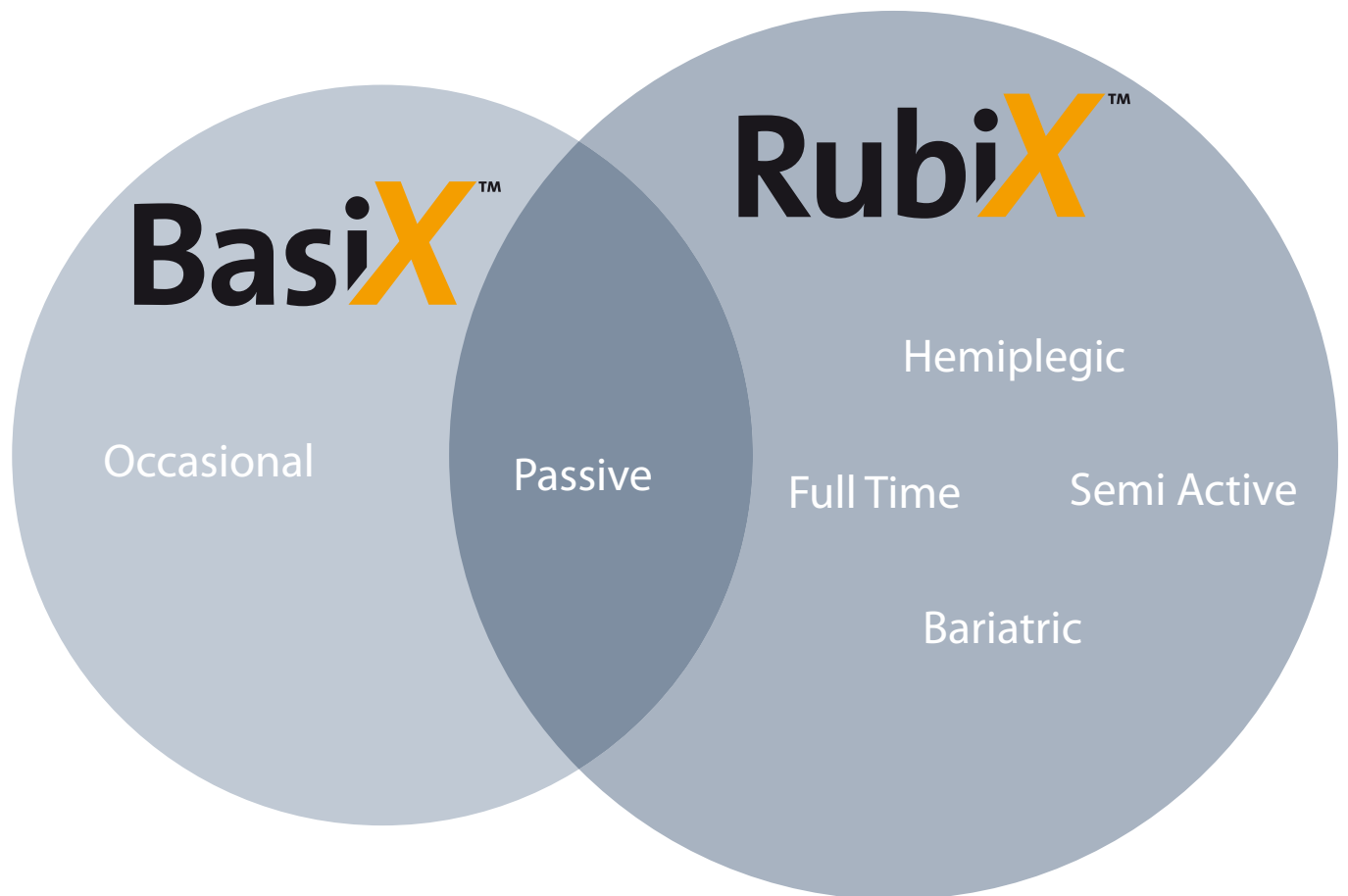

Configurability

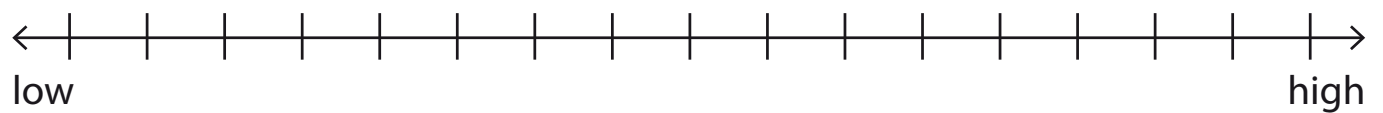

RubiX™



Full Time Use




The Breezy RubiX has been designed to ensure the optimum comfort of any full time user.

- Adjustable seat depth from 41 cm-46 cm
- Adjustable seat angle from 0° to 6°
- Adjustable seat height 42.5 cm to 50.5 cm
- Tension adjustable backrest extendable from 41 cm to 46 cm
- Length and height adjustable armrests
- Three colour options: silver, blue or red

RubiXTM

Semi Active Use

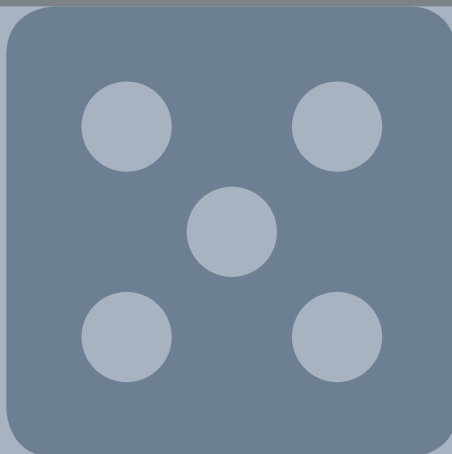

The Breezy RubiX can even be adapted to the needs of users with slightly more active lifestyles.

- The backrest can be lowered to 41 cm for a more comfortable, fluid propelling motion
- Angle adjustable seat from 0° to 6° for increased stability
- Up to 5 cm rear wheel centre of gravity adjustment for increased manoeuvrability
- Transport weight of less than 10 kg

RubiXTM

Hemiplegic Use



The Breezy RubiX can also be easily adapted to fit the needs of hemiplegic users.

- Height adjustable seat as standard from 42.5 cm. Extra-low height of 37.5 cm achievable with the short castor fork option and 6 inch castors
- Up to 5 cm rear wheel centre of gravity adjustment for increased manoeuvrability

Additional options include:

- Hemiplegic armrest
- One-arm drive
- Unilateral brake

RubiXTM

**Bariatric
Use**


**160 kg
Maximum
User Weight**



For larger, heavier users the Breezy RubiX comes in an XL version.

- Three seat width options – 52 cm, 56 cm and 60 cm
- Four seat depth options – 41 cm, 43.5 cm, 46 cm and even 48.5 cm
- Footboard and stabiliser bar for increased rigidity and support
- All key parts are interchangeable with standard BasiX and RubiX products

RubiX™

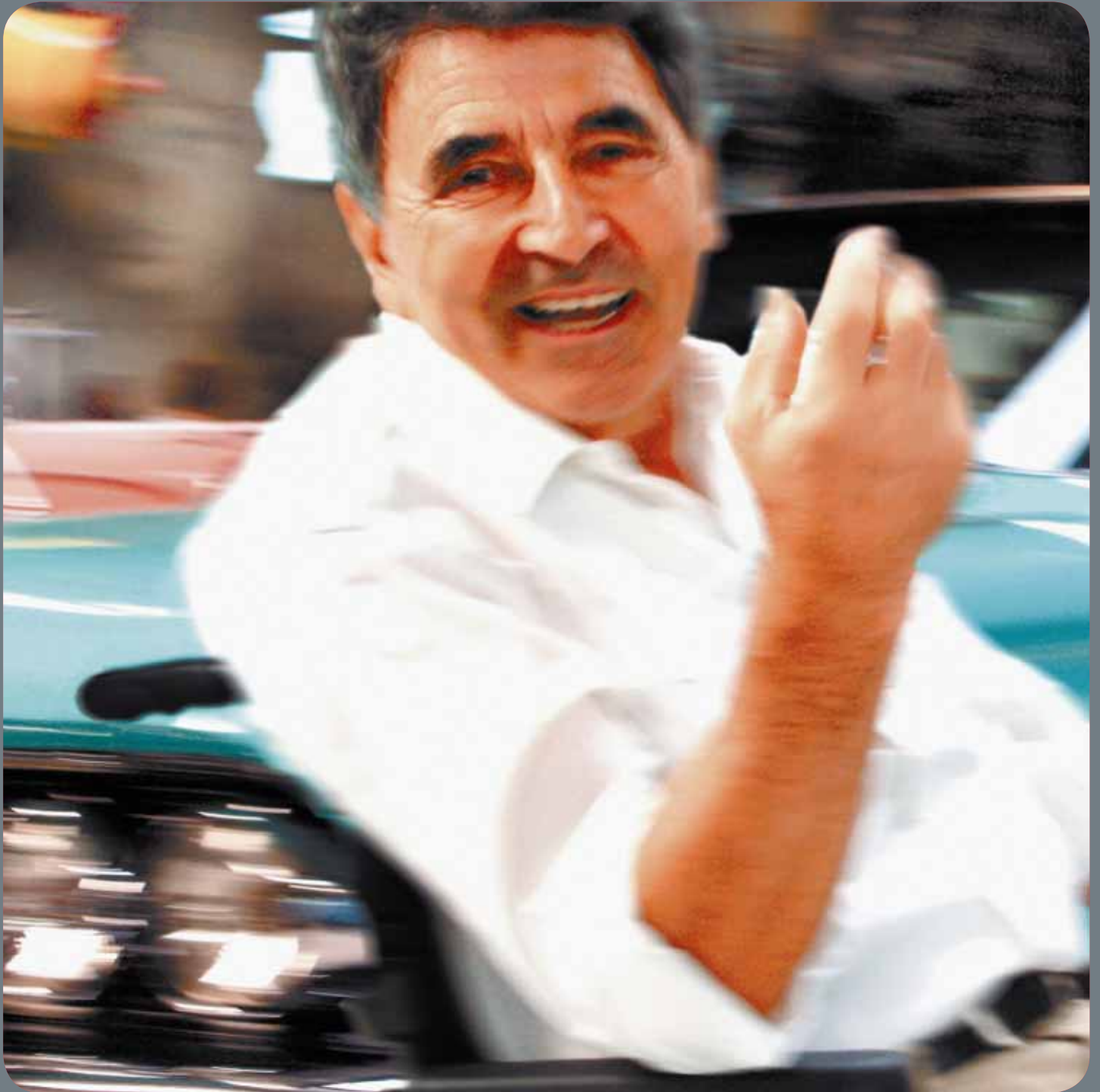

BasiX™



Passive Use



The Breezy BasiX can be quickly adapted to meet the needs of passive users

- Easy conversion to attendant push chair with optional 12 inch rear wheels and drum brakes
- Adjustable seating system that can grow with the user

Options include:

- Elevating legrest
- Amputee support
- Flip-up anti-tips

For passive users with slightly higher needs, the Breezy RubiX also has:

- Wider seat range options
- Reclining backrest
- Height adjustable push handles

BasiX[™]

Occasional Use



The comfortable and easy to use Breezy BasiX is ideal for occasional, short term or fleet use.

- Armrests, seat, backrest and footrests can be easily adjusted to meet a variety of user needs
- Tension adjustable back as standard for increased user comfort
- Optional folding back for transport and ease of storage

BasiXTM






Play the Combinations!

- Refurbishable
- Durable
- Interchangeable
- Adjustable
- Transit tested
- One chair system!
- One spares and accessories package!



Adjustable

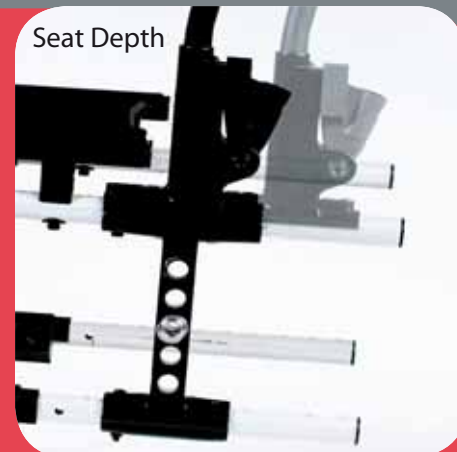

The BasiX and RubiX family is the ideal choice for those looking for a simple solution to multiple user requirements.

Both products come with the following adjustable parts as standard meaning that they can be easily modified without new parts:

- Armrest length/height
- Seat height
- Seat depth
- Seat angle
- Footplate angle
- Centre of gravity
- Backrest height
- Backrest tension



Seat Angle



Backrest Tension



Footplate Angle

Both products also use interchangeable parts so with one stock package you can service both chairs and satisfy multiple user needs quickly and easily.



Backrest Height



Seat Height

Common Spare Parts



Elevating Legrest
(without footplate) - 064260-000



Step Tube - 064931
Crutch Holder - 064970



Height Adjustable Armrest - 064475-100



Flip-up Anti-tips - 064723



Amputee Support -
064280-000



Slide On Tray - various



Tie Down Brackets - 064988

Transit Tested



Sunrise Medical takes great care in ensuring that all of our products meet the necessary testing standards and we are proud to declare that both the BasiX and the RubiX have been crash tested and fulfil the performance requirements for ISO 7176-19.

For transportation, the BasiX and RubiX should be secured using 4-part (webbing) tie down that conforms to ISO 10542 or SAE J224,9 with non-adjustable front and rear straps.

BREEZY
BasiX™
BREEZY
RubiX™

Durable and Reliable

Not only are the BasiX and RubiX flexible enough to meet most user needs, you also have the added reassurance that the product has been designed to last.



Aluminium Parts

Both products comprise of entirely aluminium components meaning that key parts such as the backrest will not wear or flex over time.



Castor Angle


No complex assemblies required. Our simple, yet durable, two bolt aluminium design can be easily adjusted and secured into place.



Three Arm Cross Bar

The additional cross bar arm leads to increased support and rigidity. Furthermore, the easy fold mechanism stays perfectly in line when unfolding and locks into place safely and securely.

Family Size Range

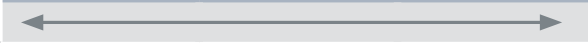
Back height	41	43.5	46
			
Seat height	37.5**	42.5	50.5
			

		Width							
		38	40.5	43	45.5	48	52	58	60
Depth*	41								
	43.5								
	46								
	48.5								

BasiX



RubiX



RubiX XL



* RubiX seat depth can be preselected to any of the listed increments or adjusted after purchase

** Achievable with 6 inch castors and short fork

Easy to Order

Not only are the BasiX and RubiX easy to use, service and adjust, we have also simplified the ordering process.

Our clear and concise order forms can be completed easily and quickly – the BasiX order form is just one page!


There are also a variety of ways you can place your order. Please see our post, e-mail or fax details below:

Sunrise Medical UK
High Street
Wollaston
Stourbridge
West Midlands, DY8 4PS
Telephone: 01384 446622 Fax: 01384 446644
E-Mail: orders@sunmed.co.uk or enquiries@sunmed.co.uk


Placing your spare parts orders could not be easier either. Thanks to our Sunparts Online service, you can place spares orders via the web and access the most up to date parts manuals.

Please visit: <http://parts.sunrisemedical.eu> for details.


BasiX™

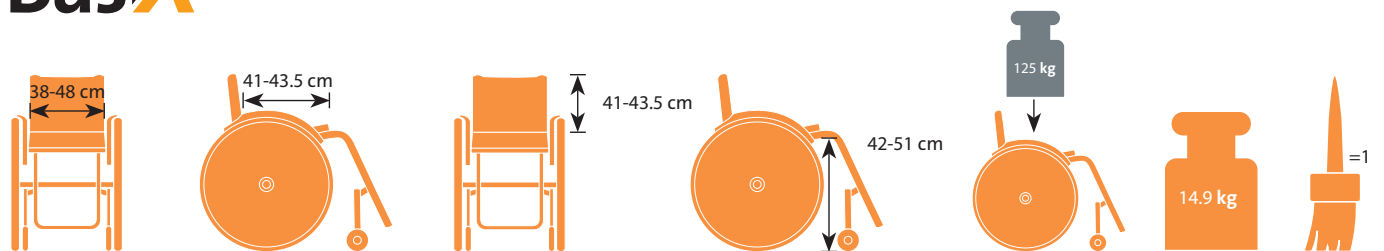

RubiX™

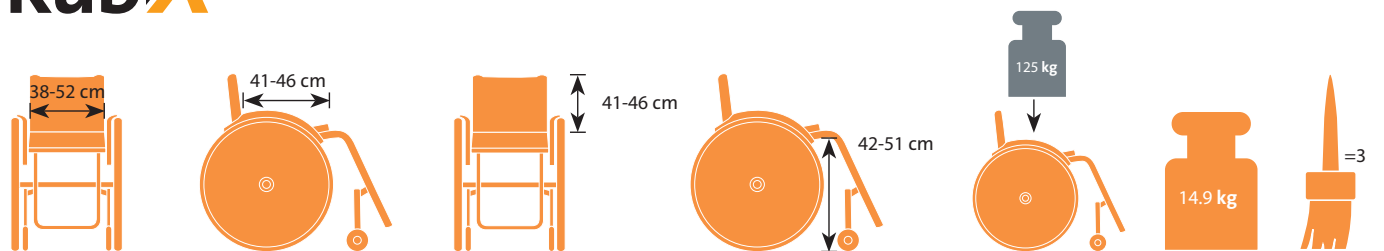

RubiX™ XL

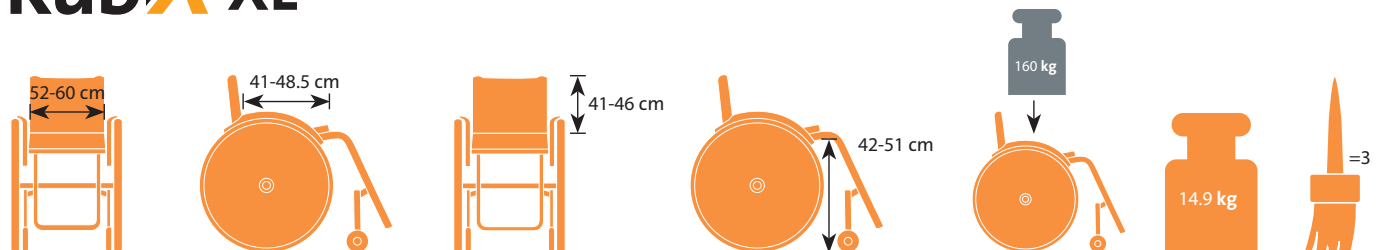



www.sunrisemedical.com