



QUICKIE
Live without limits.

XENON®

The lightest folding wheelchair in the world

XENON[®]

Minimalistic oval frame combined with an innovative, patented flat cross brace design gives the Quickie Xenon an impressive stylish and active look.



The lightest **FOLDING** wheelchair in the world! **8.4 kg*** of high performance technology.



Laura Rampini - The first paraplegic parachutist in history.

Quickie® has combined the latest ultra-light, 'hydroformed aluminium' technology, design and engineering to create Xenon - **the lightest folding wheelchair in the world.**

At just **under 8.4 kg***, and with innovative cross brace folding technology, Xenon is easy to fold, lift and transport. It has also been cleverly designed to create the kind of minimalistic style that you would expect from a rigid chair.

With high end design and innovative technology, Xenon provides a lightweight, **energy efficient driving and handling experience.**

Xenon looks and performs like a rigid chair with the benefits of a folding chair - perfect...



* 38 x 40 cm; low frame; slim style sideguards; 24" Proton wheels + aluminium handrims; 3" castors; angle adjustable back; 35 cm fixed height backrest; compact wheel lock light-weight; no pushhandles + lightweight backrest upholstery



Easy folding and lightweight lifting. Perfect for an active lifestyle.



Weighing just under **6.0 kg** without its quick release wheels, the ultra-light Xenon has revolutionised wheelchair lifting and transporting.

Incorporating the **new patented Dynamic Folding System and the Quickie® Locking System**, Xenon can be folded in seconds with just one hand.

Xenon locks automatically in the folded position with the Quickie® Locking System. To unfold simply press one lever and the Dynamic Folding System starts the unfolding process.

The cross brace has been cleverly designed to sit neatly under the seat, leaving a minimalistic, open frame that has all the style of a rigid chair.



Wheelchair shown has half folding back option

Xenon is legally protected to ensure it remains unique for you ...

Issued intellectual property rights
1080 618
001 080 733
EP 180 9532 A2

Pending intellectual property rights
PCT/ EP2008/ 011015
EP 0915 1773.0
EP2008/ 011017
EP10196521
EP10195017
DE102011 000 030.5


Laura Rampini and Filippo Landi are Quickie® Ambassadors and travel around the world with their Xenons. For more information about their amazing life and achievements please visit: www.quickiewheelchairs.eu/xenon



Experience more freedom and independence...with Xenon -
Live without limits.

XENON®



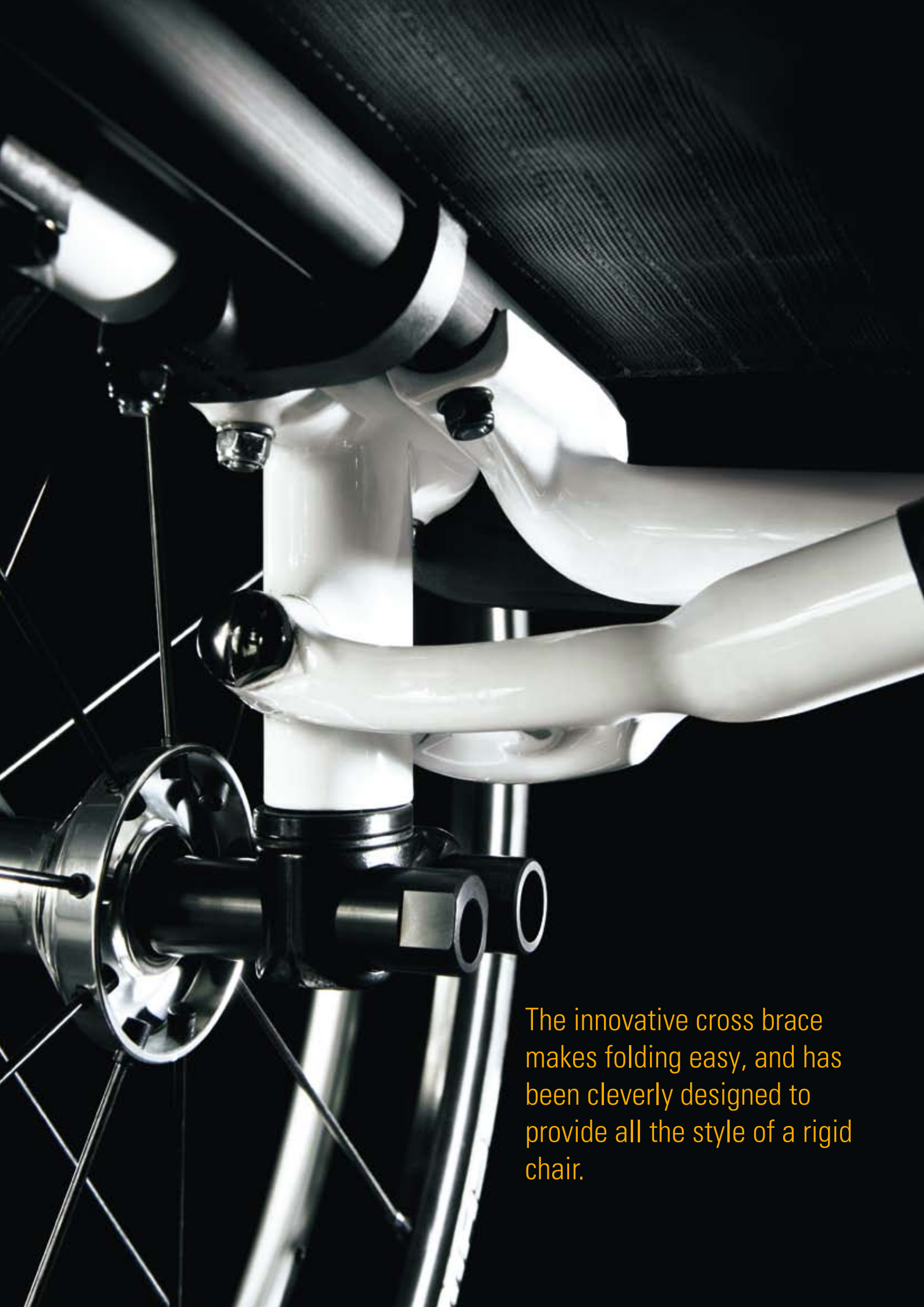

The innovative cross brace makes folding easy, and has been cleverly designed to provide all the style of a rigid chair.

Select from a wide range of options to create your own Xenon and pick from more than 20 colours to customise each component in any combination. Please visit the colour configuration tool at www.quickiewheelchairs.eu/xenon to create your personalized Xenon.



Quickie® Locking System

The innovative and patented Quickie® Locking System makes Xenon easy to handle, provides the smallest possible folding package and can also be used by people without finger function.



Lightweight back angle adjustment

Xenon offers stylish designer backrest solutions integrated into the flow of the frame. 90° fixed and angle adjustable backrests that are easy to adjust from -15° degrees to +18° in 3° increments measured to the frame.



Lightweight axle adapter

Loosen only one screw each side, turn and switch the asymmetric axle sleeves from the right to the left side and achieve a 2 cm centre of gravity change.



Standard axle adapter

Offers two centre of gravity positions, camber 0°, 2°, 4° and adjustable wheel spacing - all in one lightweight piece with quick adjustability with only two screws.




Auto-folding footrest

Provides you with the easiest / simplest folding process as it folds automatically. Lightweight and with modern looks with rounded ends, it is available with a composite and a carbon fibre plate.



Platform footrest

Available with the Xenon performance (aluminium), composite and carbon fibre footboard. It offers the lightest available footrest solution and offers full transfer space with the one side flip-up mechanism.

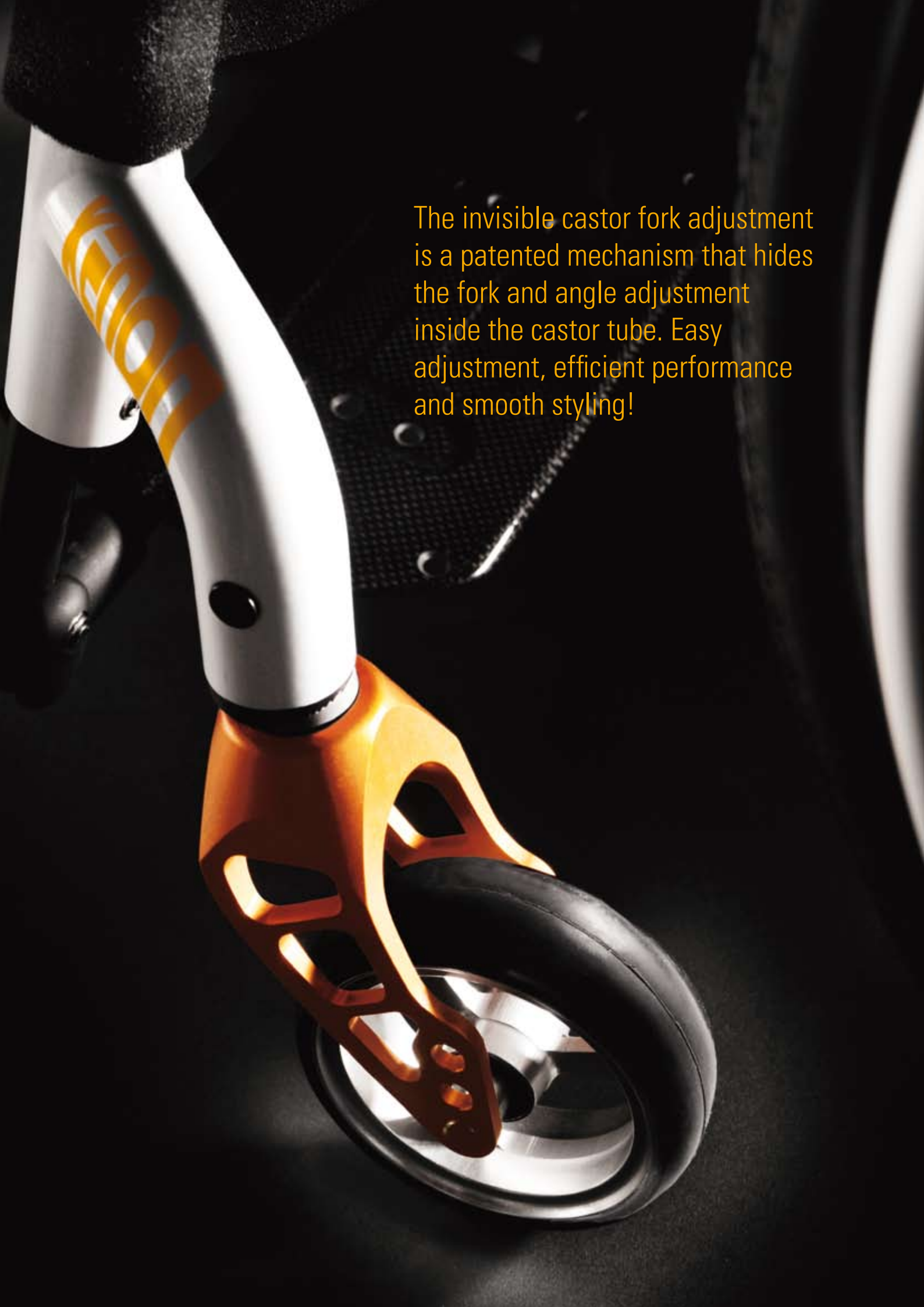A close-up photograph of a bicycle's front fork and wheel. The fork is white with yellow and white curved graphics. The front fork is orange with a lattice-like structure. The wheel has a black tire and a silver rim. The background is dark and out of focus.

The invisible castor fork adjustment is a patented mechanism that hides the fork and angle adjustment inside the castor tube. Easy adjustment, efficient performance and smooth styling!

Choose from a wide range of lightweight and comfort options to create your unique, personalised chair.



Half-folding back

Make it even easier to transport with an optional half-folding back. With this option you get the smallest folded package.



Innovative back upholstery range

Choose lightweight upholstery to reduce the weight by 130 g, or opt for the 'three piece' vented or 'Alcantara' comfort upholsteries. Breathable, adjustable, comfortable and very stylish!



Dynamic Folding System

Supports easy one hand opening after pressing one lever.



Ergonomic compact wheel locks

Minimalistic and easy to lock, the wheel locks can be screwed directly to the frame (on 0° and 2° wheel camber only). Also suitable for quadriplegics.



Proton quick release wheels

The combined weight of the wheel, handrim, axle and Schwalbe Speedrun tyre is just 2.25 kg per pair - saving over a kilogram on traditional wheels.




Suspension Fork: Frog Legs®

The Frog Legs® fork absorbs vibrations and shocks that occur constantly when riding a wheelchair. This gives the user more comfort and avoids exhaustion and long term pain issues.

XENON®

Technical data

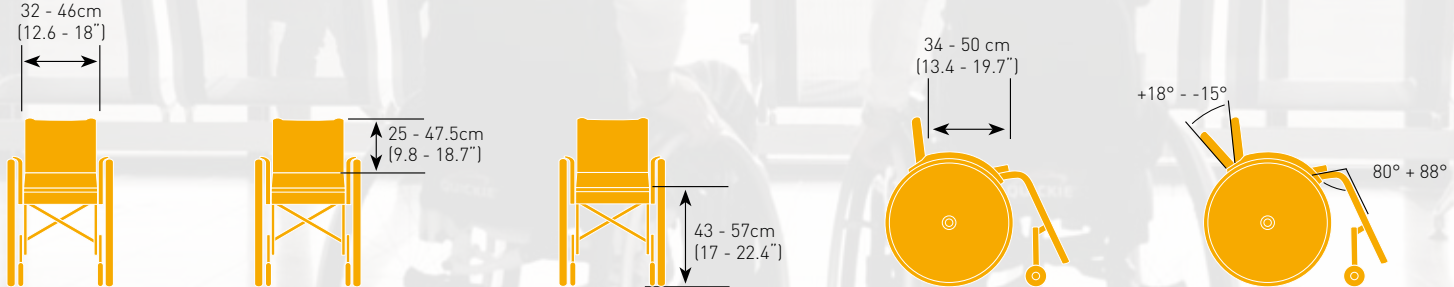

Seat width	32 - 46cm (12.6 - 18")	Seat height rear	40 - 50cm (15.7 - 19.7")
Seat depth	34 - 50cm (13.4 - 19.7")	Seat angle	0° - 12°
Back height	25 - 47.5cm (9.8 - 18.7")	Frame angle	80° & 88°
Back angle	-15° to +18°	Weight	from 8.4kg
Seat height front	43 - 57cm (17 - 22.4")	Max. user weight	125kg (19.6 stone)

To achieve the lightest possible weight with your Xenon you should choose:

- Low frame, 3" castors, 98mm x 32mm fork
- Fixed height (040227), angle adjustable backrest (040036)
- Without push handles (040038)
- Lightweight backrest upholstery (040231)
- Slim style carbon fibre sideguards (050111)
- Platform footrest (060204) with carbon footboard (060214)
- Proton wheels (080313)
- Speed Run tyres (080508)
- Aluminium handrim (080600 or 080654)
- Lightweight compact wheel locks (080803)



PB_XENON_EN_GB_REV1

Sunrise Medical Limited
 High Street, Wollaston
 Stourbridge
 West Midlands
 DY8 4PS
 Tel: +44 (0) 1384 44 66 66
 Fax: +44 (0) 1384 44 66 44
www.SunriseMedical.com