

# ARGON<sup>2</sup>

## THE LIGHTWEIGHT ADJUSTABLE RIGID CHAIR

The logo for QUICKIE, featuring a stylized orange swoosh above the word "QUICKIE" in a bold, sans-serif font. Below it is the tagline "Live without limits." in a smaller, lighter font.

QUICKIE<sup>®</sup>  
Live without limits.



The logo for MMS MEDICAL, with "MMS" in a large, bold, blue font and "MEDICAL" in a smaller, blue font below it.

**MMS**  
MEDICAL

# ARGON<sup>2</sup>

New design and more adjustments for users with higher expectations.



# ARGON<sup>2</sup> – A STRONG AND RELIABLE DESIGN FOR AN ACTIVE LIFESTYLE.

MAX USER WEIGHT  
**140**  
KG \*

Compact frame:  
knee angle 88° or 80°



## // KEEP ACTIVE – A WHEELCHAIR YOU CAN TRUST.

Personalise your Argon<sup>2</sup> with two different knee angles and frame insets to meet your requirements of a compact wheelchair. The open frame design combined with the low weight of the chair ensures easy transfer of the chair across the body into a vehicle.

The new Argon<sup>2</sup> is a strong, energy efficient and rigid wheelchair with a maximum user weight of 140 kg. It's approved for handbike usage and is crash tested according to ISO 7176-19 for transportation in a vehicle.

\* HD option must be selected

The new Argon<sup>2</sup> is the ideal combination of a rigid performing and multi-adjustable wheelchair. These functional characteristics are built in the modern looking open frame design. The Argon<sup>2</sup> is the perfect companion for your active lifestyle.

## > Handbike approved - strong and reliable

Argon<sup>2</sup> is approved for handbike attachment. You know you have a chair that's strong and reliable when it's handbike approved.



# OUR FINEST TECHNOLOGY.

Cascaded down from our flagship product the Quickie Helium.

Hydroformed castor housing. Stylish, easier to adjust and lighter for better performance.



New - Carbotecture forks, lighter and stronger

New backrest angle – -31° to +15°, lighter and easier to adjust

Infinitely adjustable centre of gravity. Slides & clamps on the frame



Camber options of 3° / 6° / 9°

Rear seat height & seat angle adjustments (1cm increments)

## // ARGON<sup>2</sup> – USING HELIUM TECHNOLOGY FOR IMPROVED ENERGY EFFICIENCY

Our engineers are always looking for innovative technologies to improve wheelchair performance and design.

The outstanding design of the castor link tube (thanks to tube hydroforming) helps to achieve complex tube geometries - for a more modern design, additional strength and better rolling efficiency.

The Argon<sup>2</sup> forks are made of Carbotecture material – a special lightweight and extremely strong carbon fibre material from the aeronautic industry. Extremely impact resistant and break proof, weight efficient (tensile strength/density) and 50% lighter than aluminium, making it the right material for our forks.

# ARGON<sup>2</sup> – FINE TUNE YOUR PERFORMANCE.

The wide range of adjustments allows you to fine tune your performance. It's the ideal wheelchair for users who are developing their driving capabilities over time and require more active chair settings. This ensures the best driving performance now and in the future.



- > **Hydroformed castor link tube**  
Modern design with additional strength and lighter

## // THE CHAIR WHICH CAN BE CONSTANTLY ADAPTED TO YOU

The Argon<sup>2</sup> offers full adjustability without compromising on performance. The rear seat height is adjustable in 1cm steps so the seat angle can be constantly adapted to changing needs. With the infinitely adjustable centre of gravity you can quickly change your driving settings to a more active or a more stable performance.

The backrest angle is also adjustable, as well as the backrest height and the tension of the back upholstery - so you'll always find the seating position that best suits you.

- > **The innovative Carbotecture fork**  
Made out of stiff high-tech material and extremely lightweight at half the weight of the aluminium forks. Available as a style package in three different colours.



Front seat height adjustments

Options and features of the new Argon<sup>2</sup> – Configure your wheelchair to your individual needs. For further information on product options please check the Argon<sup>2</sup> order form.



**// CARBON SIDEGUARD**

Very individual and lightweight sideguard, with or without fender.



**// CARBON PLATFORM FOOTPLATE**

Lightweight and stylish footplate with angle and depth adjustment.



**// ANTI-TIP FLIP-UP**

Easily operated by the user, this flip-up design ensures easier transfer of the chair across the body into a vehicle.



**// ONE ARM WHEEL LOCK**

Easy to use one hand operated wheel lock.



**// TOOL HEIGHT ADJUSTABLE ARMREST**


Provides easy transfer and good arm support with a low transport weight.



**// PADDED TUBULAR ARMREST**

Due to the swing-away possibility this armrest gives you support when needed and does not hinder during propulsion.

ACTIVE, MODERN, RELIABLE.  
PERFECT FOR YOU.

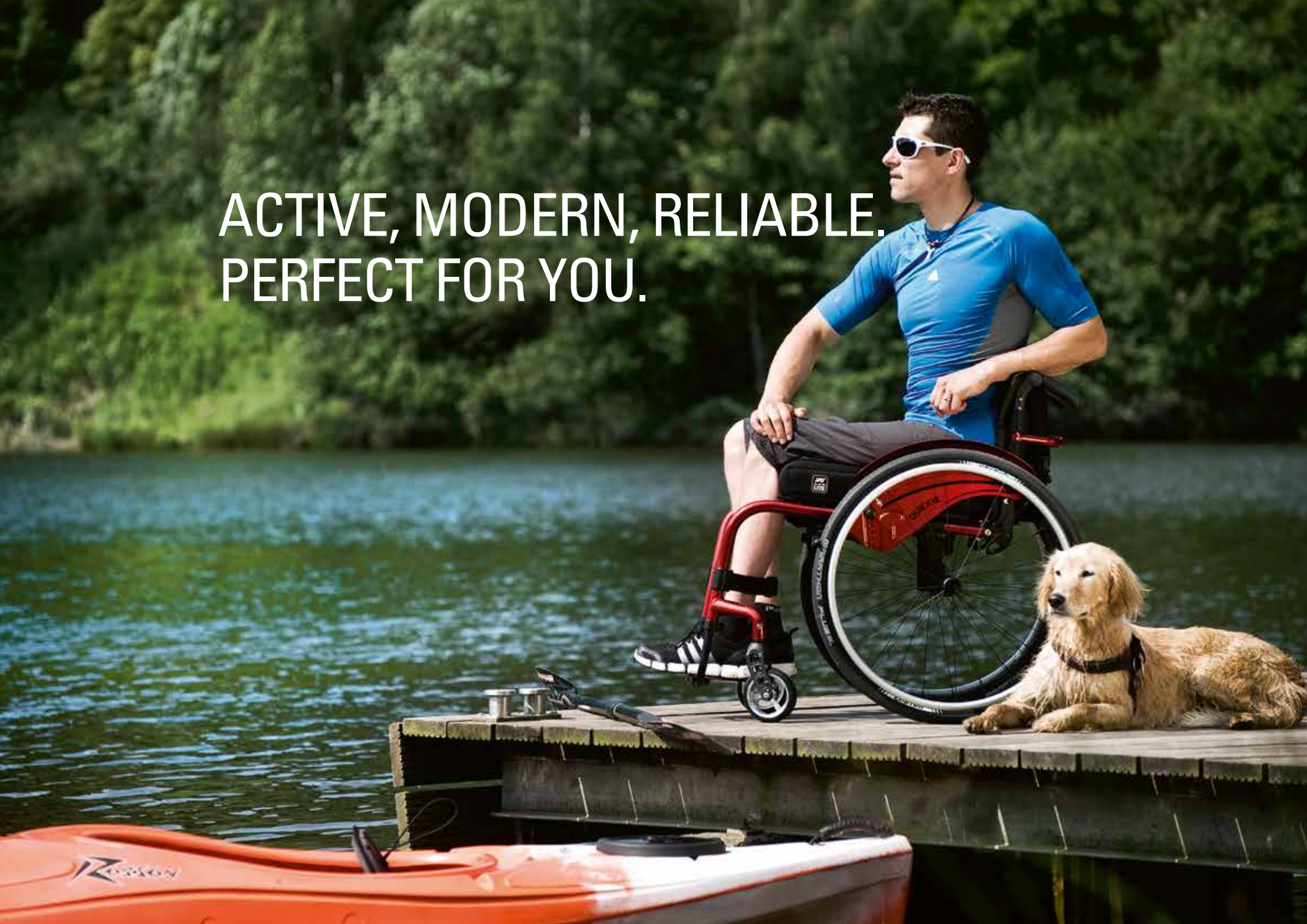

# TECHNICAL DATA

## ARGON<sup>2</sup> // Rigid wheelchair

Max. user weight: 125 kg / 140 kg

▸ Proton wheels or handbike adapter: 100 kg

▸ With lightweight wheels: 125 kg

Seat width: 300 – 500 mm

Seat depth: 300 – 500 mm

Seat height (front): 450 – 520 mm

Seat height (rear): 340 – 520 mm

Back height: 250 – 450 mm

Backrest angle: -31° to + 15°

COG: 32 to 104 mm

Camber: 0°/ 3°/ 6°/ 9°

Frame angle: 88° or 80°

Inset: 0 mm or 30 mm

Chair weight: from 9 kg

Colours: 32 colours

Crash tested: Yes

Find out more today!

Visit [www.SunriseMedical.co.uk](http://www.SunriseMedical.co.uk) for further information on product options, or contact your local dealer or representative.



### Get SOCIAL with QUICKIE

Keep up-to-date with the latest news on the Quickie Argon<sup>2</sup> and other Quickie products through our Facebook page:

[www.facebook.com/QuickieWheelchairs](http://www.facebook.com/QuickieWheelchairs)



MMS Medical, 51 East Gate Drive, Little Island, Cork

Tel: 021 4618000 Fax: 021 4618099 E: [info@mmsmedical.ie](mailto:info@mmsmedical.ie)

MMS Medical, Unit 81 Western Parkway Business Park, Lower Ballymount Road, Dublin 2

Tel: 01 4500582 Fax: 01 4565889 E: [sales@mmsmedical.ie](mailto:sales@mmsmedical.ie)

[www.mmsmedical.ie](http://www.mmsmedical.ie)

