



QUICKIE
Live without limits.

XENON²
SUPERIOR FOLDING AND
DRIVING PERFORMANCE

MMS
MEDICAL

XENON² – THE STYLISH, HIGH-PERFORMANCE SERIES

Whether you're looking for the very best in driving dynamics or the most compact folding dimensions, there's a Xenon² just for you.

Will it be the clean, streamlined design of the Xenon² FF? Alternatively how about the reinforced rigidity of the Xenon² Hybrid? Or the practical portability of the Xenon² SA? Whichever one of the three feature-packed models you choose, the sophisticated elegance of the Xenon² will always make you the centre of attention.



” Music is her passion. Her voice is expressive: Monika Kuszyńska represented Poland in the Eurovision Song Contest in Vienna performing her song “In the Name Of Love”. In Poland she is already very successful and is also sharing her message with people outside the country.



XENON² FF THE LIGHTEST

Lightest folding wheelchair of the world.
Open fixed front frame.
Cool, clean and streamlined design.



XENON² HYBRID THE STRONGEST

No compromises in terms of rigidity.
Reinforced fixed front frame.
Maximum user weight 140 kg.



XENON² SA THE ALLROUNDER

Easy standing transfers in and out of
your wheelchair.
Reinforced frame with swing away footrests.
Most compact folded dimensions for easy
portability by car.





THE LIGHTEST FOLDING CHAIR IN THE WORLD...

How do you achieve a folding wheelchair that weighs as little as 8.8 kg? By using premium frame materials like aerospace-grade aluminium and an innovative process that's unique only to Quickie – ShapeLoc Technology.

By heat-treating the aluminium, we're able to increase its strength, which allows for the frame's tube walls to be thinner, lighter and stronger than a conventional folding wheelchair. Just feel it...





...WITH PEERLESS FOLDING PERFORMANCE

At the core of the Xenon² is a unique cross-brace that fits so neatly under the seat it's barely noticeable – giving a very minimalistic, open-frame feel that you'd usually associate with a rigid wheelchair.

But the similarities don't end with the appearance – the completely re-designed axle stem, combined with the new cross-brace, provides the stiff, stable driving performance feel of a rigid frame, but with all the portability of a folding one.

→ **Easy to use folding mechanism**

The new aluminium latches on the cross-brace will also keep your wheels in parallel when folding, so your Xenon² remains stable, no matter the situation.






STYLE IN MOTION

Personalise your Xenon² and stand out from the crowd with a versatile palette of 32 frame colours. There's also 5 distinctive anodised styles to choose from for the castor fork and wheel. You can even accentuate your colour scheme further by matching or contrasting the upholstery trim colour of the EXO backrest upholstery.

→ **Personalize your Xenon²**
with a selection of 32
eye-catching colours

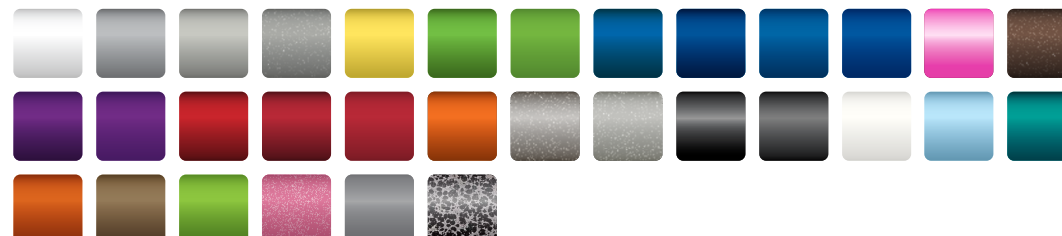

” **Monika Kuszyńska:** I like fashionable outfits. The attractive design of the Xenon² compliments my creative style perfectly.



→ **Add some colour to your back upholstery**
This very comfortable and stylish Exo back upholstery can be chosen with 5 different colour binding tapes.



→ **Choose your castor and fork own style**
The anodized front castors and the forks are selectable in silver, blue, red orange or black for an unique look



” **Monika Kuszyńska:** *My Xenon² is very like my music, it feels like a part of me. It's like a companion when I travel all around the world.*

ADJUSTABLE TO SUIT YOU

Exactly as it should be – designed to be adaptable to your changing needs, the axle stem on the Xenon² provides easy height adjustment.

You want to change your chair driving characteristic from passive to more active or the opposite way? Just change the centre of gravity setting by simple repositioning of the axle stem. And don't forget to adjust the back angle for your ergonomic posture.

1 Centre of Gravity

Enables a driving characteristic from very active to passive.

2 Backrest Angle

The right back angle adjustment is important for your ergonomic posture.

3 Rear Seat Height

Exactly as it should be - adjust the height to your personal need.





XENON[®] FF
THE LIGHTEST



XENON[®] HYBRID
THE STRONGEST



XENON[®] SA
THE ALLROUNDER

XENON²

Ultra lightweight wheelchairs
with outstanding folding and
driving performance



Options and features of the Xenon² series – Configure your wheelchair to your individual needs. For further information on product options please check the Xenon² and Xenon² SA order form.



// FOLD-DOWN BACKREST

This half-folding backrest folds to the front and achieves a small pack size of the wheelchair for transports.



// SURGE[®] LT HANDRIMS

Its Gripton rubber strip on top of the ergonomic oval shaped rim offers maximum traction.



// ANGLE ADJUSTABLE BACK

Stylish designer backrest solution integrated into the follow of the frame. Backrest angle from 75° to 103° in 4° increments.



// DESK SIDEGUARD

This flip-up and removable sideguard can be ordered with fixed or height adjustable arm pads.



// PADDED TUBULAR ARMREST

Due to the swing away possibility this armrest gives you the support when needed and do not hinder during propulsion.



// CARBON SIDEGUARD

Very individual and lightweight sideguard with its unique look.



// AUTO-FOLDING FOOTBOARD

This footplate is angle and depth adjustable and available with carbon or composit. Available for fixed front frame versions.



// TIE-DOWN BRACKETS FOR TRANSPORTATION

With this option the wheelchair is approved for transportation in a vehicle (crash tested according to ISO 7176-19).

TECHNICAL DATA

	Xenon ² FF	Xenon ² Hybrid	Xenon ² SA
Max. user weight:	110 kg	140 kg	125 kg
Seat width:	32 – 46 cm	32 – 50 cm	32 – 50 cm
Seat depth:	34 - 50 cm	34 - 50 cm	34 - 50 cm
Seat height front:	43 – 55 cm	43 – 55 cm	43 – 56 cm
Seat height rear:	37 – 50 cm	37 – 50 cm	37 – 50 cm
Backrest height:	25 – 47,5 cm	25 – 47,5 cm	25 – 47,5 cm
Backrest angle:	75° to 103° in 4° increments	75° to 103° in 4° increments	75° to 103° in 4° increments
Centre of gravity (COG):	4 – 13 cm	4 – 13 cm	4 – 13 cm
Camber:	0° / 2° / 4°	0° / 2° / 4°	0° / 2° / 4°
Frame Version:	Open frame with fixed front	Reinforced frame with fixed front	Reinforced frame with swing-away front
Chair weight:	from 8,8 kg	from 9,5 kg	from 10,3 kg
Frame colours:	32	32	32
Crash tested:	yes	yes	yes

For further information on the full specification, options and accessories, please refer to the order form.

All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.

If you are visually impaired, this document can be viewed in PDF format at www.SunriseMedical.eu or alternatively is available on request in large text.



MMS Medical, 51 East Gate Drive, Little island, Cork.
Tel : 021 4618000 Fax: 021 4618099 E: info@mmsmedical.ie

MMS Medical, Unit 81 Western Parkway Business Park, Lower
Ballymount Road, Dublin 2
Tel: 01 4500582 F: 01 4565889 E: sales@mmsmedical.ie



www.mmsmedical.ie