

F5 Corpus VS



Get vertical.

Stability, Comfort, Performance.

The F5 Corpus VS is a durable, powerful and stable front wheel drive power wheelchair combined with the best seating system on the market – the Corpus seat. The seating system is designed to effectively support the clinical, functional and lifestyle needs of active users.

The latest version of the Vertical Stander (VS) offers all the power seat functions of the Corpus seat plus an innovative new way of standing, allowing you to choose from two positions, Sit-to-stand or Lay-to-stand.

Giving you greater independence and freedom, the F5 Corpus VS features fully independent suspension and automatic, self-adjusting support wheels, providing a new level of comfort and stability, especially when you stand and drive, indoors or outdoors!

A range of alternative driving controls for non-joystick drivers is available, as well as separate controls for assistants/caregivers.

For transport purposes, the upper part of the backrest can be removed to reduce the overall size of the wheelchair.

permobil

THE POWER OF MOBILITY

electrical functions





Stand and drive

Choose between Sit-to-stand or Lay-to-stand, or have both! Stand and drive mode is also made easier at any standing angle.

Comfort

All the benefits of the ergonomically designed Corpus seating system are now available on the Corpus VS.

Stability

Powered self-adjusting support wheels for stability while standing and safe negotiation of obstacles while standing and driving.

Control


The latest driving and seat positioning controls. Full colour LCD display with memory functions, speedometer and clock.

Colours



sunburst orange radiant red power pink ocean blue titan blue galactic green phantom grey midnight black

Upholstery

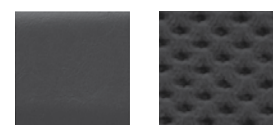

leatherette black fabric

Specifications

Length	1145 mm	Legrest angle – adjustable	90°–180°
Width	655–790 mm	Battery capacity	2 x 73 AH
Height	960–1170 mm	Range	25–35 km*
Min. transport length	980 mm	Max. speed	12 km/h
Min. transport height	825 mm	Turning diameter (ISO 7176-5)	1525 mm
Seat height	450 mm	Turning in corridor	1140 mm
Seat height with electric seat lift	450–800 mm	Obstacle capability	70–100 mm
Weight incl. batteries	196 kg	Suspension	Yes
Weight of batteries	2x23 kg	Electronics R-net	120A
Max. user weight "Sit-to-stand" (max 80°)	136 kg	Electric seat elevator	Yes
Max user weight "Lay-to-Stand" (max 80°)	110 kg	Electric tilt adjustment	50°
Seat width	420–570 mm (by 50 mm)	Electric legrest adjustment	Yes
Seat depth	370–570 mm (by 50 mm)	Electric backrest adjustment	Yes
Backrest width	360/410/460/510 mm	Electric stand-up function	Yes
Backrest height	595/620/645/670 mm	Warranty	2 years
Distance between armrests	380–480/480–580 mm	ISO 7176-19/ISO 10542-3	Yes – crash tested
Backrest angle – adjustable	180°	EN12182/EN12184 CE	Yes
Armrest height	215–330 mm		

These specifications are dependent on the configuration of the wheelchair. For further information not mentioned in this sheet, please refer to the user manual which can be downloaded from the Permobil web site, or contact your local Permobil office. *The range may be reduced if the wheelchair is used frequently on slopes, rough ground or to climb kerbs often.