



Blizzard



Multi Adjustable Rear Wheel Drive Power Chair

The **Easicare Blizzard** is a successful combination of user lead design and high quality build. The overall result is a highly adjustable power chair with excellent all-round performance for use outdoors and indoors, together with user comfort and durability.

Excellent Performance

The **Easicare Blizzard** utilises the Dynamics Electronic Control System for control and manoeuvrability with two motor rear wheel drive for performance. With a range of 16km / 10 miles and a top speed of 8kmh / 5mph it provides enhanced mobility without compromise.

Mechanical Adjustability

Inherent adjustability is included in the design specification resulting in easy and quick to adjust settings to accommodate different user needs.



Fold-up and swing-away footrests



Kerb climber (Optional)



Head rest (Optional)



Height adjustable armrests



Easy to reach charging socket



Positioning belt as standard



Anti tip wheels fitted for safety



13" Alloy rear wheels



Fully folding back



Adjustable tilt positioning




Additional tilt adjust



Dynamics Electronic Control System as standard



Adjustable backrest as standard

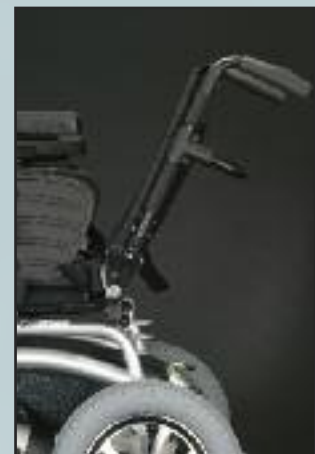



AVAILABLE OPTIONS:

- Transport headrest
- Tray
- Attendant control
- Swing away controller
- Kerb climber
- Elevated leg rest
- 40Ah battery upgrade

FEATURES:

- Maximum Speed: 8kmph / 5 mph
- Travel Range: 20km / 12 miles (36Ah);
25km / 15.5 miles (40Ah);
- Capacity: 136Kg / 300lbs
- Durable Structure
- Back Angle Adjustable
- Height Adjustable Armrests
- Left / Right Side Exchangeable Controller
- Length Adjustable Controller Assembly
- Adjustable seat width / depth / height

SPECIFICATIONS:

- Overall Length: 1040mm / 41"
- Overall Width: 580mm / 23"
- Overall Height: 920mm / 36"
- Wheels - Front: 200mm / 8"
- Wheels - Rear: 320mm / 13"
- Weight with Batteries: 74Kg / 163lbs
- Ground Clearance: 100mm / 4"
- Grade Climbable: 8 degrees
- Kerb Climbing: 30mm / 1 1/4"
- Turning Radius: 920mm / 36"
- Brake: Electro-Mechanical
- Seat Width: 355 - 560mm / 14" - 22"
- Seat Depth: 405 - 460mm / 16" - 18"
- Backrest Height: 405 - 510mm / 16" - 20"
- Backrest Angle: 100 - 120 degrees
- Drive Train: 2-Motor Rear-Wheel Drive
- Battery Weight: 24Kg / 52 lbs
- Motor Size: 420 W 4600 rpm
- Battery: (2) 12V. 36Ah
- Charger: 5A Off Board
- Electronics: Dynamic 50A
- Seat Type: Padded, breathable ballistic nylon back rest



MMS Medical, 51 Eastgate Drive, Little Island, Cork
Tel : 021 4618000 Fax : 021 4618099 Email : info@mmsmedical.ie

MMS Medical, Unit 81 Western Parkway Business Park, Lower
Ballymount Road, Dublin 12.
Tel : 01 4500582 Fax : 4565889 Email: info@mmsmedical.ie

www.mmsmedical.ie

For further information on the full specification of options and accessories of the Easicare Blizzard Power Chair, please refer to the order form.
All information is subject to change without notification.
Please consult your dealer with any queries you may have. E&OE.