

# WHILL Model C

“The next-generation powerchair that is redefining accessibility.”

It all starts with great design. Your powerchair should enable you to be spontaneous and free but to do it all in style.

Use the simple iPhone app and Bluetooth technology to drive your WHILL Model C remotely, bring it over when you need it and store it out of the way when you don't. Lock and unlock your WHILL using the app and monitor all functions including battery life. You can even store and charge your phone onboard.

One of the most striking features of the WHILL Model C is the omni-wheel. Like no other wheel, it allows you to turn on the spot with the tightest turning circle on the market. There are no traditional castors so you have complete 360° control. It's tough too. Go up and down curbs with ease and travel on gravel, grass and uneven surfaces with confidence.

From the ergonomic controller to the electromagnetic brakes, the WHILL Model C puts you firmly in command. With fully responsive systems and four speed settings, your WHILL is agile and reactive.

The WHILL Model C is designed to complement your lifestyle so it goes wherever you go. It disassembles in just 10 seconds into three lightweight pieces that easily fit into a standard car boot.

- Unique omni-wheel – for on the spot turning
- Bluetooth technology – for remote operation
- Intuitive controls – for safe and agile movement
- Ultra-portable – 10-second disassembly into three compact parts
- Stylish design – sleek contours and choice of six colours
- Fully adjustable – to fit your lifestyle

[www.mmsmedical.ie](http://www.mmsmedical.ie)

The WHILL Model C is the next generation powerchair that's shifting the boundaries of mobility. From indoors to outdoors it combines unmatched precision and manoeuvrability with the strength and adaptability to take on the most unpredictable terrain. In one chair you can navigate the smallest room with ease and move straight outdoors to tackle curbs, pathways and trails with confidence.

Full details at [www.tgamobility.co.uk/whill](http://www.tgamobility.co.uk/whill)

**MMS**  
MEDICAL



# WHILL Model C




Maximum Speed (up to)	Maximum Weight	Range (up to)	Total Weight
4mph	16 stone	10 miles	115kgs


## The unique omni-wheel system



The innovative alternative to traditional caster wheels, the large omni-wheels have 24 revolving sections to allow on the spot rotation within a very tight 30" (76cm) turning circle.

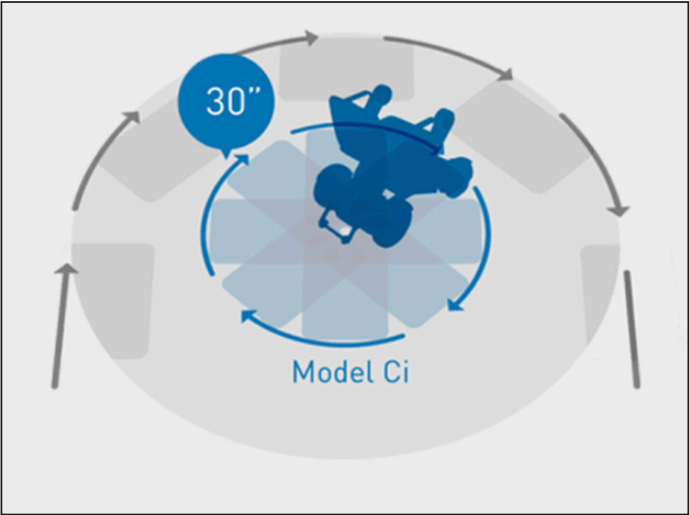


## The right fit

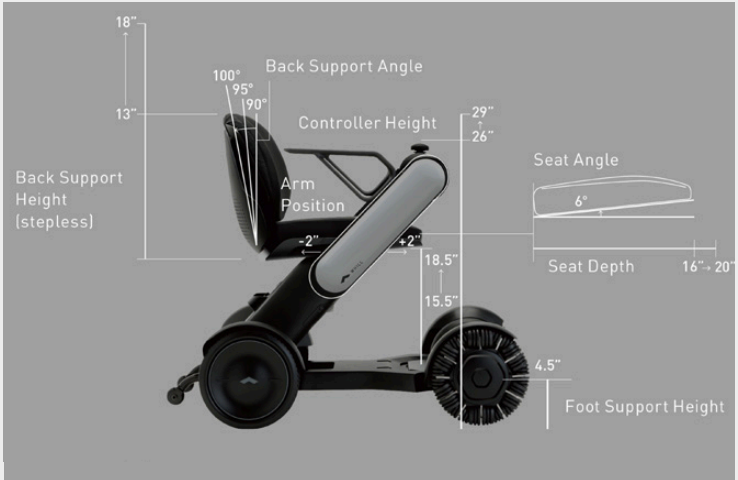

Model C can be easily adjusted from the seat depth and height to the controller position (right or left drive), to fit you and create a comfortable driving experience.

## Advanced ground clearance

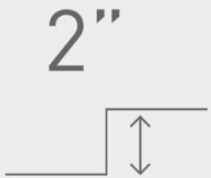

Effortlessly climb over obstacles up to 2" (5cm) in height with its two powerful motors and large front omni-wheel.




Travel unhindered over rough terrain with a 2.25" (5.7cm) ground clearance.



**Ergonomic controls**

Comfortable controls for effortless steering and safe manoeuvrability. Including four speed settings, a safety horn and battery life indicator.

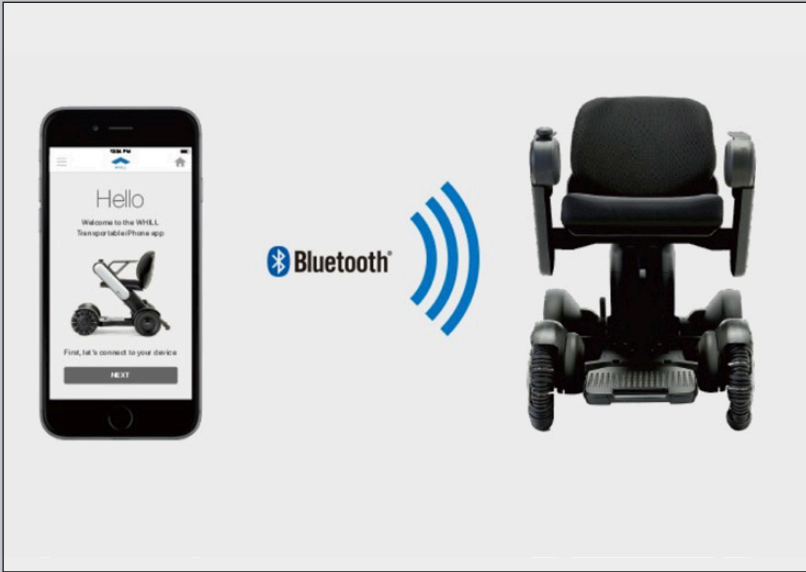

**Technology at your fingertips**

Featuring remote control driving and secure lock and unlock function via an easy to use Apple iOS app. Monitor performance and battery life in real time.



**Lightweight, long-life battery**

With a 10-mile range and only five hours charge time.



**Seat options**

WHILL Model C comes with the choice of either a 16" or 18" seat.

**Feel safe and secure**



Anti-sway feature helps keep the device driving in a straight line on side slopes.



Electromagnetic brakes provide smooth and safe stopping, even on inclines.



Footplate folds up and arm rests easily rotate out of the way for easy access and barrier free use at tables.



Voice notifications alert you when Model C encounters operational errors.





### Ultra-portable design

It takes just 10 seconds to disassemble the WHILL Model C into just three compact pieces. The entire chair fits easily into a standard sized boot.



### Multi award winner

#### The Red Dot Design Award

An internationally recognised product design and communication design prize awarded by the Design Zentrum Nordrhein Westfalen in Essen, Germany.



**reddot award**  
best of the best  
2019

- iF Design Award
- CES 2018 Innovation Award
- International Design Excellence Award
- Japan Good Design Award

### Find your style

There is a choice of six individual metallic colours to personalise the WHILL C.



### Accessories

WHILL C has a great selection of optional accessories to match your personal style and to enhance your active lifestyle.



Underseat bag



Lap belt



Smart key



Walking stick holder



Tablet holder



Side bag storage and arm

**MMS Medical**  
51 East Gate Drive,  
Little island  
Cork

Tel : 021 4618000  
Fax: 021 4618099  
Email: [info@mmsmdical.ie](mailto:info@mmsmdical.ie)  
[www.mmsmedical.ie](http://www.mmsmedical.ie)