




Proudly made
and manufactured
in Canada!

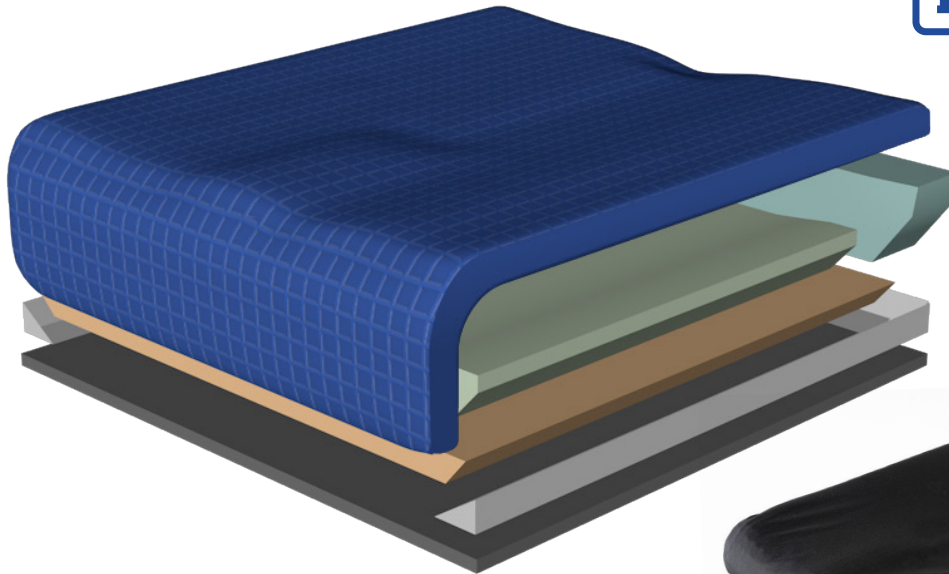

Product Specifications

CONSTRUCTION

Closed cell foam base, dual high-density foam core with a firm sacral shelf and cascading Pillow Top gel topper

COVERS

Two outer cover options: Breathable Spacetex fabric with incontinent protection OR Elastatex soft touch outer incontinent cover.

Inner cover: Elastatex breathable, wipe-able, antibacterial, fluid proof cover.

STANDARD SIZES

40x40, 40x45, 40x50,
45x40, 45x45, 45x50,
50x40, 50x45, 50x50,
55x40, 55x45, 55x50,
55x55 cm

CUSTOM SIZES

All dimensions available

PRODUCT CODE

GHLP"WxD" eg. GHLP4545

WEIGHT CAPACITY

Recommended maximum 115kg

*Can be upgraded to 160kg weight capacity upon request

The Geo-Matrix™ Hybrid Low Profile cushion has the same benefits as the standard Hybrid, however, with a 2.5cm height reduction. Making it an ideal solution for foot propulsion & obtaining the ideal seat to floor height.

The Geo-Matrix™ Hybrid Low-Profile cushion is often used to provide comfort for users who require gentle positioning and contouring with excellent pressure distribution and vibration control.

The cascading gel top layer increases comfort and reduces pressure behind the knees. The Geo-Matrix™ Hybrid Low-Profile is designed for moderate postural support.

Blake Medical esteemed distributor:



T: 021 461800
E: info@mmsmedical.ie
www.mmsmedical.ie




WASHING INSTRUCTIONS

Covers: 60°C in washing machine. Hang to air dry. Do not bleach.

Geo-Matrix™ gel and foam: Wipe down gel using any cleaner that does not contain bleach and allow to air dry.

For more specific cleaning instructions please visit our website or user's manual.

All of Blake Medical products are free from phthalates and latex. This includes all medical fabrics, Geo-Matrix™ gel, moulded polyurethane, standard polyurethane and specialty foam. All cushions meet EN-1021 part 1 & 2 requirements for fire retardancy.

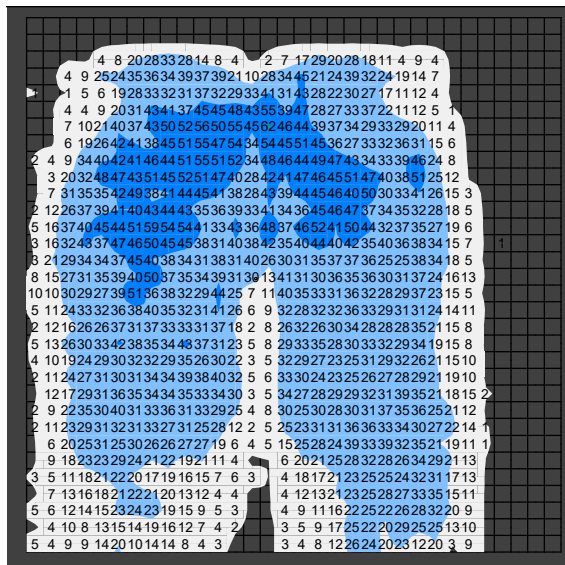

Client: Male 70kg

ADDITIONAL BENEFITS

- Cushion remains 2-3 degrees cooler under seated body temperature to assist with microclimate management.
- Minimizes shear through column segmentation.
- Geo-Matrix™ Gel does not change structure in climates of extreme heat nor extreme cold.
- Geo-Matrix™ gel is recyclable through Blake Medical.

Advantages of Geo-Matrix™ Gel

