

Data Sheet

300H

The 300H is our lightweight, easy-to-handle evacuation chair, designed for safe and comfortable use with passengers up to 182kg. Its double-handed and anti-slip grips, combined with an optimised hammock, provide stability and comfort during evacuation. The unique friction belt system ensures smooth, controlled movement on stairs, making the 300H our most popular choice for single-person operation.

Specification

Height	1100mm (folded)	1350mm (unfolded)
Width	520mm (folded)	520mm (unfolded)
Depth	210mm (folded)	810mm (unfolded)
Weight (unladen)	9.5kg	-
Payload Capacity	182kg	-
Stair Angle	28° - 40°	-

MMS
MEDICAL



EVAC+CHAIR®

Emergency Evacuation For Everyone

Features

- 1. Double handed grip**
Enhances user operation.
- 2. Contoured headrest and integrated safety strap**
Provides passenger support and comfort.
- 3. Optimised hammock**
Greater passenger comfort.
- 4. Sewn-in safety belt**
Improves safety and durability.
- 5. Simplified graphical instructions and user guide**
Easy to follow step-by-step instructions.
- 6. Optimised kickstand, with repositioned rear wheels**
Improves stability and enhances operator handling.
- 7. Upgraded castor wheels**
Strong 50mm wheels improving manoeuvrability.
- 8. Unique friction belt system**
Self-braking track controls smoothness and speed of descent.
- 9. Anti-slip grip**
Improves stability and operation during passenger transfer or when deploying the kickstand.

Approvals/Accreditations

- + TUV NORD
- + CE & UKCA Mark
- + Registered Class 1 Medical Device
- + ISO 13485: Medical Device Quality Management
- + ISO 9001: Quality Management System

Benefits

- + Easy to operate
- + Single person operation
- + Lightweight
- + No heavy lifting involved
- + 23-point quality check throughout manufacture

Includes

- + Photoluminescent signage
- + Wall mounted hooks
- + Dust Cover
- + User guide



300H

Accessories

Discover our range of must-have accessories – engineered to boost safety, comfort and confidence. Whether upgrading your Evac+Chair 300H or adding extra support for passengers, make sure you're fully equipped to protect every life.

Comfy Seat

Provides passengers with additional comfort and support during training and evacuation. Designed and manufactured from memory foam.



Secure Steel Cabinet

Keep your Evac+Chair safe and secure with a lockable steel cabinet, featuring an acrylic window and customisable left or right opening (subject to order quantity).




Dust Cover

Protect your Evac+Chair from dust, pollutants and damage with our flame-retardant Dust Cover.



Evac+Chair Stand

For quick and easy access, safely and securely store next to the stairs or in the designated area.

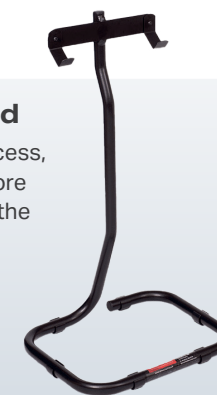

Reception Notice

Inform your visitors that there is an Evac+Chair available onsite, offering reassurance that in an emergency, they can be safely evacuated.



Photoluminescent Sign

To be positioned above your Evac+Chair, our Photoluminescent Sign will help locate your Evac+Chair in poor visibility conditions.




Anti-theft Device

Protecting your Evac+Chair will reduce the chance of theft, damage or misuse. Our anti-theft device acts as a deterrent and can be easily installed in a matter of minutes.



Transfer Board

Ensures a safe and efficient transfer from a seat to an Evac+Chair for mobility impaired individuals. The Transfer Board serves as a temporary bridge, facilitating a smooth and seamless transition.

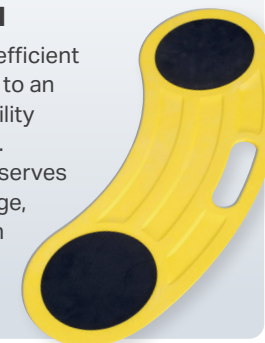

T +44 (0)121 706 6744
E enquiries@evacchair.co.uk
W evacchair.co.uk

T :021 4618000
E : info@mmsmedical.ie
W : www.mmsmedical.ie