


beka 
hospitec

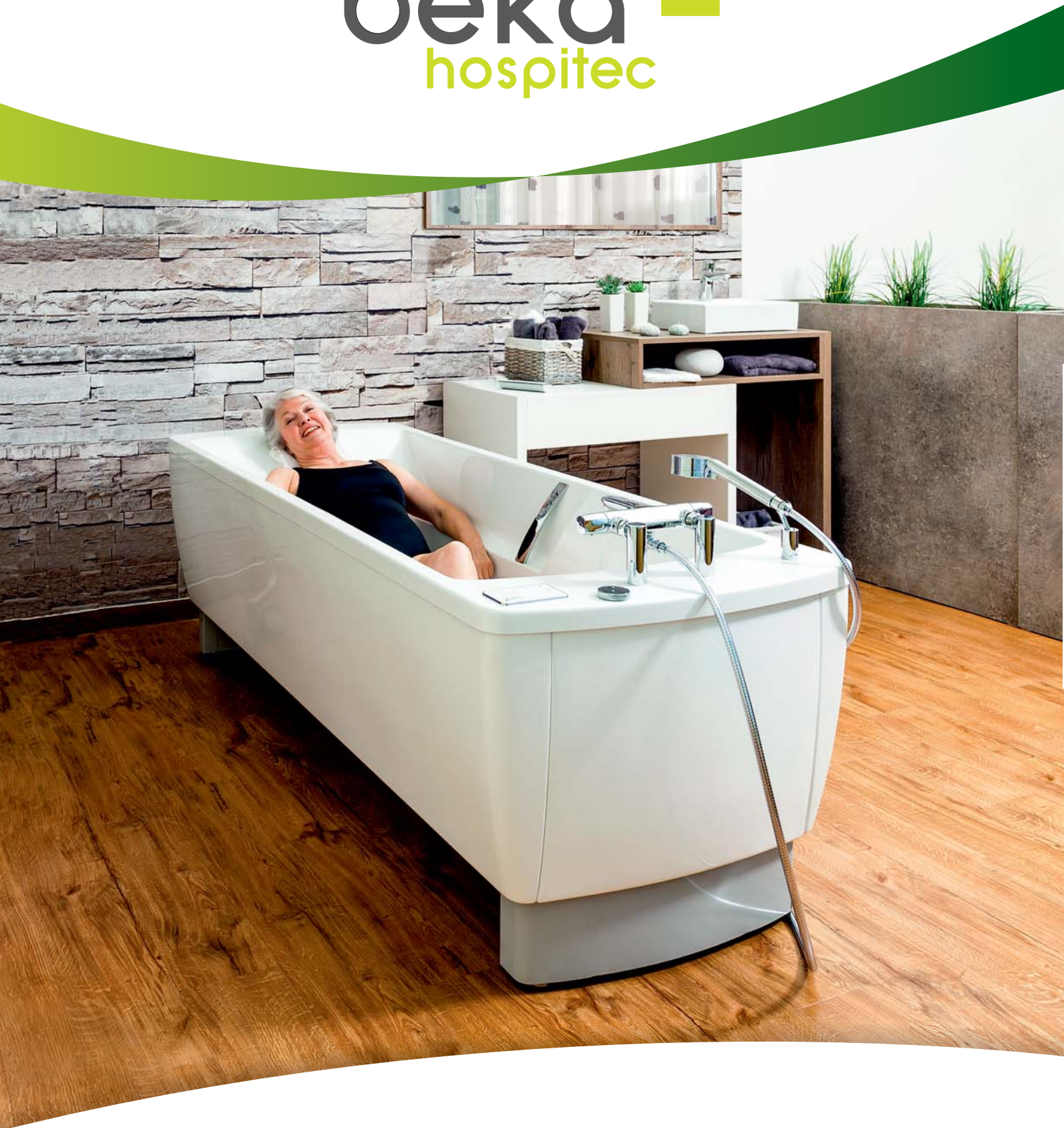

AVERO Comfort

Wellness and care tub

MMS
MEDICAL

AVERO Comfort lifting tub

Elegant design for your nursing bathroom ...



Designed by Uwe Spannagel™

Optional change of the viewing direction ...

The clear and soft shape of this tub reminds of modern, contemporary private bathrooms and perfectly fits in the existing wellness and bathing concept of your bathroom.

A special feature of the AVERO Comfort tub are its radii which are reflected both on the inside and on the outside lining. Without compromising the beauty of the tub or the aesthetics of your nursing bathroom, this unique tub design enables the tub to be rotated and placed in the room in such a way that the resident can look freely around the room while bathing. This greatly enhances the feeling of wellbeing of the bather.





Suitable for all care hoists in the market ...



The AVERO Comfort lifting tub stands solidly in the room without fixation and can be freely accessed so that even in constricted spaces and in the narrowest of spaces, all modern care hoists available in the market can be used without any problems, regardless of the manufacturer. The size of the comfortable interior enables the patient to lie down comfortably and to bathe with or without hoist.

The specially adapted back shape of the tub body with its prominent arm rests makes for a harmonised interior tub form which allows for ergonomic lying comfort to the patient. A thermostatic mixer of the Hans-Grohe brand mounted on the tub rim strengthens the homely character of the lifting tub.

AVERO Comfort lifting tub

Clear contours and optimum functionality ...



Designed by Uwe Spannagel™

Designer tub with homely character ...



The modern, elegant and timeless tub shape fits unobtrusively in your new nursing bathroom and at the same time gives your residents, your staff and visitors the impression of a domestic surrounding.

Premium, chromed brass grips with distinctive design and additional insertions for a tub shortener embedded in the tub rim not only give the tub a distinctive look but also provide excellent hold points for the patient.



Clear contours and optimum functionality ...



Options for the AVERO Comfort


Bathing and wellness of the highest level ...

A full bath has beneficial effects on the body and the mind. According to the latest findings in dementia research, sensory impulses, in particular, have shown a therapeutic effect on demented residents ...



Coloured light therapy system ...

The coloured light therapy system of the AVERO Comfort creates a special and optically stimulating atmosphere of wellbeing for your resident. Two LED lights built into either side of the bath tub provide a relaxing and therapeutic bathing experience without blinding the resident unnecessarily. The coloured light therapy system has a freely selectable colour option and an automatic sequence programme for colour sequences. We also offer a superior coloured light system for the AVERO Comfort.

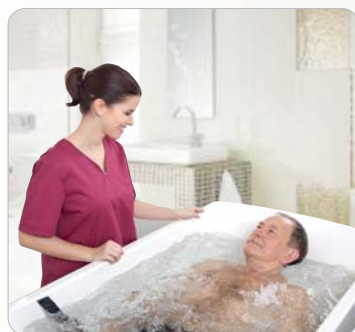

Air spa system ...

An integrated air spa system with its cleansing and vitalising abilities provides a comfortable massaging effect which is highly beneficial to the resident. The slight tingling on the skin is relaxing and soothing, especially in case of restless patients, and further promotes the blood circulation of the skin. Contamination of the integrated jets in the tub bottom and the foot area is not possible due to a special process, and there is therefore no need for an integrated disinfection system to disinfect the air spa system.



Music sound system ...

With two sound pressure generators integrated in the area of the upper body part, the vibrations of the music are transmitted directly into the water and to the resident's body. This creates a tactile sound experience which has a positive influence on the bather, as shown by the latest research. A USB port allows to play the music chosen by your resident. It is operated via wireless remote control.

When it comes to colour, you can select from a variety of RAL colours for our AVERO Comfort. Of course, it is also possible to select different colours for the inside and the outside of the tub ...



Aromatherapy system ...

An aromatherapy system integrated in the tub body allows for different aromas to create a wonderful experience for the senses of the resident in conjunction with the air heated by the air spa system and the bath water. It goes without saying that the operation is simple and the system is maintenance-free. This option is available with the air spa system only.



Change of the viewing direction ...

Without compromising the beauty of the tub or the aesthetics of your nursing bathroom, this unique tub design enables the tub to be rotated and placed in the room in such a way that the resident can look freely around the room while bathing. This greatly enhances the feeling of wellbeing of the bather.



Automatic water stop ...

Take advantage of the option to fill the tub with an automatic filling stop to be able to pay full attention to the resident whilst the tub is being filled. You can set the filling level as you wish so that you can select the right amount of water adapted to the resident and the application. This saves time and energy costs and allows for the nursing staff to fully concentrate on the resident. Operation is by means of a keypad.

The optimum solution for your nursing staff and your residents ...



AVERO Comfort combines wellness and nursing perfectly and supports your nursing staff optimally in the daily care routine ...

Benefits for your nursing staff:

- ❑ Easy operation by means of a keypad and symbols
- ❑ Continuous height adjustment for optimum working height
- ❑ Non-porous surface for easy cleaning and optimum hygiene
- ❑ The optional automatic water stop saves time and costs and allows for the nursing staff to fully concentrate on the resident



The AVERO Comfort lifting tub offers your resident and your nursing staff maximum safety and comfort at all times ...

Benefits for your resident:

- ❑ Ergonomic interior tub body for comfortable and safe bathing with and without hoist
- ❑ Arm rests with grips for more comfort and safety
- ❑ Optional tub shortener for a secure hold when bathing
- ❑ Domestic appearance
- ❑ Scalding protection for hand shower and tub spout
- ❑ Numerous options for relaxing wellness bathing

AVERO Comfort lifting tub

Wellness and care tub ...





Benefits for investors and utilities management:

- ❑ Modern design for a modern nursing bathroom
- ❑ Low spare parts and maintenance costs
- ❑ Proven technology, not dependent on a factory customer service
- ❑ Use of branded components from reputed manufacturers, e.g., motors, valves, control systems
- ❑ Hans-Grohe branded thermostat
- ❑ Very good price-performance ratio
- ❑ Free technical support over the phone by BEKA Hospitec

Technical data of the AVERO Comfort

Overall dimension, incl. operation console	2.100 mm x 800 mm
Inside tub dimension (L x W x D)	1.770 mm x 575 mm x 400 mm
Lifting range	680 - 1.080 mm
Lifting capacity	4000 N per drive
Frequency	50 Hz
Power consumption	max. 4,5 ampere
Power supply	230 Volt
Patient weight (max.)	250 kg
SWL (max.)	400 kg

Accessories / Options

Change of viewing direction in the bathroom
Automatic water stop
Tub shortener
Air spa system
Coloured light therapy system
Superior coloured light therapy system
Music sound system
Aromatherapy system
Special colours as per RAL



MMS Medical Ltd.
51 East Gate Drive Little
Island,
Cork.

Tel: 021 4618000
Fax: 021 4618099
E: info@mmsmedical.ie

MMS Medical Ltd.
Unit 81, Western Parkway Business Park,
Lower Ballymount Road,
Dublin 12.

Tel: 01 4500582
Fax: 01 4565889
www.mmsmedical.ie