

beka hospitec



AVERO Phoenix

Wellness and care tub in an exclusive design

MMS
MEDICAL

AVERO Phoenix lifting tub

Clear design and high-quality material ...



Designed by Uwe Spannagel™

Tub body with lowered foot area ...



The AVERO Phoenix tub is freely accessible from underneath on three sides so that it can be easily used with bath lifts even in restricted spaces in very narrow rooms. The tub stands firmly on the floor with no fixation and can be optimally adjusted to your nursing bathroom. The tub body and foot area of the AVERO Phoenix tub allow for the use of all modern bath lifts and hoists available in the market.

The **lowered leg area** and the **specially adjusted back profile** of the tub body, with distinctive arm supports, form a balanced ergonomic inner shape of the tub which allows the resident to lie down in the tub in a relaxed and ergonomic manner. This proves to be of great advantage especially when using a seat lift, as it provides additional support to the resident inside the tub. The seating position is also equipped with a hygienic opening section which facilitates easy access for intimate cleaning.



Safety and comfort while bathing ...

The operating and control unit of the AVERO Phoenix lifting tub impresses with a distinctive elaborate shape and icons which eliminates the risk of faulty operation and facilitates work for staff. Special safety systems for the tub spout and hand shower protect the resident from scalding and inadvertent activation of the disinfection facility as long as the resident is lying in the tub.

High-quality chrome-plated grab handles in an unique design and additional slots integrated in the tub rim for the tub shortener give the tub not only a distinctive look but also excellent support and a feeling of safety to the resident who can enjoy the bath in complete relaxation.



AVERO Phoenix lifting tub

The optimum solution for your nursing staff and your residents ...



Designed by Uwe Spannagel™

Bathing experience at the highest level ...



The AVERO Phoenix tub has been developed on the basis of the latest findings and in collaboration with nursing staff and occupational therapists. Selected high-quality materials, a clear design and above all the ergonomic tub interior body give the resident an opportunity to enjoy a bathing experience at the highest level.

The bathroom is given a completely new look with an appealing design. A simple and elaborate operating and safety concept allows for easy and simple operation by staff and at the same time provides complete safety for the resident and patient.



Benefits for your nursing staff and your residents ...


- ❑ Simple operation via a membrane keypad containing icons
- ❑ Safety modules prevent faulty operation
- ❑ Non-porous surface for simple cleaning and optimum hygiene
- ❑ Works in an ergonomic manner, height adjustment of the tub to the working height required
- ❑ Ergonomic tub interior body allows comfortable and safe bathing of the resident with or without hoist
- ❑ Arm supports with grab handles for added comfort and safety
- ❑ Integrally formed seat with hygiene opening section for better intimate cleaning
- ❑ More than sufficient space for the resident's shoulders
- ❑ Double scalding protection for the hand shower and tub spout
- ❑ Suitable for all hoists



Options for the AVERO Phoenix

Wellness and therapy ...

A full bath has beneficial effects on the body and the mind. According to the latest findings in dementia research, sensory impulses, in particular, have shown a therapeutic effect on demented residents ...

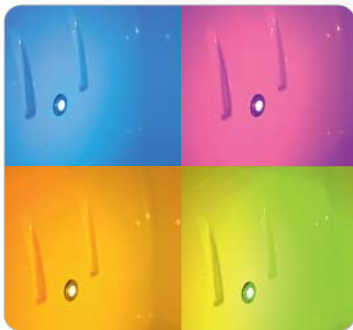

Coloured light therapy system ...

The coloured light therapy system of the AVERO Phoenix creates a special and optically stimulating atmosphere of well-being for your resident. By means of four LED lights integrated in the foot area, the bathing water can be immersed in the entire colour spectrum without heat radiation and without blinding the resident unnecessarily. The coloured light therapy system has a freely selectable colour option and an automatic sequence programme for colour sequences.



Air spa system ...

An integrated air spa system with its cleansing and vitalising abilities provides a comfortable massaging effect which is highly beneficial to the resident. The slight tingling on the skin is relaxing and soothing, especially in case of restless patients, and further promotes the blood circulation of the skin. Contamination of the integrated jets in the tub bottom and the foot area is not possible due to a special process, and there is therefore no need for an integrated disinfection system to disinfect the air spa system.



Music sound system ...

With two sound pressure generators integrated in the area of the upper body part, the vibrations of the music are transmitted directly into the water and to the resident's body. This creates a tactile sound experience which has a positive influence on the bather, as shown by the latest research. A USB port allows to play the music chosen by your resident. It is operated via wireless remote control.

When it comes to colour, you can, of course, select from a variety of RAL colours for our AVERO Phoenix and optimally adjust the colour of your tub to your furnishing concept ...



Aromatherapy system ...

An aromatherapy system integrated in the tub body allows for different aromas to create a wonderful experience for the senses of the resident in conjunction with the air heated by the air spa system and the bath water. It goes without saying that the operation is simple and the system is maintenance-free. This option is available with the air spa system only.



Automatic filling of the tub ...

Enjoy the opportunity of filling the tub with an automatic water stop whilst being able to fully concentrate on the resident in the process. The AVERO Phoenix offers two different filling levels to choose from to be able to select the right amount of water depending on the resident and the application. This saves time and energy and allows to devote the fullest of attention to the resident.



Disinfection facility ...

The disinfection facility installed offers maximum protection against cross-infections and allows the use of suitable, authorised, in-house manufactured surface disinfectant concentrates. The disinfection hand shower and storage box are secured behind the control area out of reach so that the resident cannot make contact inadvertently while bathing.

Technical data of the AVERO Phoenix

Overall dimensions AVERO Phoenix 190	2.284 mm x 890 mm
Tub interior body AVERO Phoenix 190 (L x W x D)	1.780 mm x 575/720 mm x 430/490 mm
Tub external body AVERO Phoenix 190 (L x W x D)	1.934 mm x 800 mm x 520 mm
Height adjustment	680 mm - 1.080 mm
Lifting capacity	8.000 N
Power consumption	max. 4,5 ampere
Filling capacity of the AVERO Phoenix 190	approx. 220 litres
Max. patient weight	220 kg
SWL	300 kg

Accessories / Options

Automatic water stop
Tub shortener
Disinfection facility with separate hand shower
Air spa system
Coloured light therapy system
Music sound system
Aromatherapy system
Special colour as per RAL



MMS Medical Ltd.
51 East Gate Drive
Little Island,
Cork.

Tel: 021 4618000
Fax: 021 4618099
E: info@mmsmedical.ie

MMS Medical Ltd.
Unit 81, Western Parkway Business Park,
Lower Ballymount Road,
Dublin 12.

Tel: 01 4500582
Fax: 01 4565889
www.mmsmedical.ie