

beka 
hospitec



AVERO Premium

The ideal wellness and care tub

MMS
MEDICAL

AVERO Premium lifting tub

Combines ergonomics and modern design ...



Designed by Uwe Spannagel™

Safety and comfort while bathing ...



The simple and elaborate operating and safety concept of the AVERO Premium lifting allows easy and simple operation of the tub and thus prevents faulty operations. This makes work for staff significantly easier. An integrated branded thermostat, with scalding protection for the tub spout and hand shower, protects the resident against scalding and provides the nursing staff with the necessary safety in their daily routine.

High-quality chrome-plated grab handles in a unique design and additional slots integrated in the tub rim for the tub shortener give the tub not only a distinctive look but also support and a feeling of safety to the resident who can enjoy the bath in complete relaxation.



Designed by Uwe Spannagel™

Tub body with lowered foot area ...



The AVERO Premium tub is freely accessible from underneath on three sides so that it can be easily used with bath lifts even in restricted spaces in very narrow rooms. The tub stands firmly on the floor with no fixation and can be optimally adjusted to your nursing bathroom. The tub body and foot area of the AVERO Premium tub allow for the use of all modern bath lifts and hoists available in the market.

The **lowered leg area** and the **specially adjusted back profile** of the tub body, with distinctive arm supports, form a balanced ergonomic inner shape of the tub which allows the resident to lie down in the tub in a relaxed and ergonomic manner. This proves to be of great advantage especially when using a seat lift, as it provides additional support to the resident inside the tub. The seating position is also equipped with a hygienic opening section which facilitates easy access for intimate cleaning.

AVERO Premium lifting tub

The optimum solution for your nursing staff and your residents ...



Bathing experience at the highest level ...



The AVERO Premium tub has been developed on the basis of the latest findings and in collaboration with nursing staff and occupational therapists. Selected high-quality materials, a clear design and above all the ergonomic tub interior body give the resident an opportunity to enjoy a bathing experience at the highest level.

The modern tub body together with the elegant lifting unit presents an attractive image. A simple and elaborate operating and safety concept allows for easy and simple operation by staff, and at the same time provides complete safety to the resident and patient.



Benefits for your nursing staff and your residents ...


- ❑ Simple operation
- ❑ Safety modules prevent faulty operation
- ❑ Non-porous surface for simple cleaning and optimum hygiene
- ❑ Works in an ergonomic manner, height adjustment of the tub to the working height required
- ❑ Ergonomic tub interior body allows comfortable and safe bathing of the resident with or without hoist
- ❑ Arm supports with grab handles for added comfort and safety
- ❑ Integrally formed seat with hygiene opening section for better intimate cleaning
- ❑ Scalding protection for the hand shower and tub spout
- ❑ More than sufficient space for the resident's shoulders
- ❑ Suitable for all hoists



Options for the AVERO Premium

Wellness and therapy ...

A full bath has beneficial effects on the body and the mind. According to the latest findings in dementia research, sensory impulses, in particular, have shown a therapeutic effect on demented residents ...

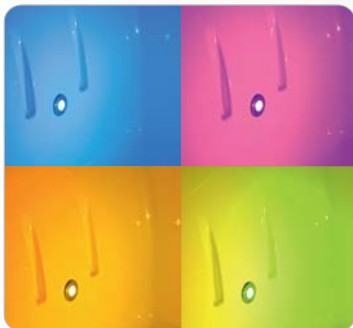

Coloured light therapy system ...

The coloured light therapy system of the AVERO Premium creates a special and optically stimulating atmosphere of well-being for your resident. By means of two LED lights integrated in the foot area, the bathing water can be immersed in the entire colour spectrum without heat radiation and without blinding the resident unnecessarily. The coloured light therapy system has a freely selectable colour option and an automatic sequence programme for colour sequences.



Air spa system ...

An integrated air spa system with its cleansing and vitalising abilities provides a comfortable massaging effect which is highly beneficial to the resident. The slight tingling on the skin is relaxing and soothing, especially in case of restless patients, and further promotes the blood circulation of the skin. Contamination of the integrated jets in the tub bottom and the foot area is not possible due to a special process, and there is therefore no need for an integrated disinfection system to disinfect the air spa system.



Music sound system ...

With two sound pressure generators integrated in the area of the upper body part, the vibrations of the music are transmitted directly into the water and to the resident's body. This creates a tactile sound experience which has a positive influence on the bather, as shown by the latest research. A USB port allows to play the music chosen by your resident. It is operated via wireless remote control.

When it comes to colour, you can, of course, select from a variety of RAL colours for our AVERO Premium lifting tub and optimally adjust the colour of your tub to your furnishing concept ...



Aromatherapy system ...

An aromatherapy system integrated in the tub body allows for different aromas to create a wonderful experience for the senses of the resident in conjunction with the air heated by the air spa system and the bath water. It goes without saying that the operation is simple and the system is maintenance-free. This option is available with the air spa system only.



Automatic filling of the tub ...

Enjoy the opportunity of filling the tub with an automatic water stop whilst being able to fully concentrate on the resident in the process. The filling level can be freely selected and is then set once. After the water stop is activated, you can always press a button and continue filling the tub. This saves time and energy and allows to devote the fullest of attention to the resident.



Disinfection facility ...

The disinfection facility installed offers maximum protection against cross-infections and allows the use of suitable, authorised, in-house manufactured surface disinfectant concentrates. The disinfection hand shower and storage box are secured behind the control area out of reach so that the resident cannot make contact inadvertently while bathing.

Technical data of the AVERO Premium

Overall dimensions AVERO Premium 190	2.230 mm x 820 mm
Tub interior body of the AVERO Premium 190 (L x W x D)	1.780 mm x 575/720 mm x 430/490 mm
Tub exterior body of the AVERO Premium 190 (L x W x D)	1.934 mm x 800 mm x 520 mm
Height adjustment	662 mm - 1.062 mm
Lifting capacity	10.000 N
Power consumption	max. 4,5 ampere
Max. patient weight	220 kg
SWL	300 kg

Accessories / Options

Choice of tap fittings
Automatic water stop
Tub shortener
Disinfection facility
Air spa system
Coloured light therapy system
Music sound system
Aromatherapy system
Special colour as per RAL



MMS Medical Ltd.
51 East Gate Drive
Little Island,
Cork.

Tel: 021 4618000
Fax: 021 4618099
E: info@mmsmedical.ie

MMS Medical Ltd.
Unit 81, Western Parkway Business Park,
Lower Ballymount Road,
Dublin 12.

Tel: 01 4500582
Fax: 01 4565889
www.mmsmedical.ie