

## PRODUCT DATASHEET

### Digital Scale for Clinical Use, Rectangle Model

Brand	: CURAS®
Category	: Nursing
Range	: Medical Scales
Product Group	: Scale & Weight
REF	: 05901



### PRODUCT DESCRIPTION

The Digital Scale for Clinical Use is specially designed for Curas Hygiene Pulp™. The rectangular Model, 20x16.5X1.9cm, is with tempered glass surface. It can weigh up to 5kg. With lighted LCD, Auto Zero & Off function. Comes with 2xAAA Batteries.

### INTENDED USE

The digital scale is specially designed for weighing and measuring medical pulp product and its content. Also, suitable for other medical items.

### FEATURES & BENEFITS

- Specifically designed for pulp urine bottles (REF 08050 - 52, 08056 - 60). For measuring and converting of urine content from weight to volume.
- Default unit, in grams, and can be converted to other measurement (i.e. lb, oz, ml).
- Able to weight at minimum 5g and maximum 5kg.
- Tempered glass surface that is easy to maintain, kept clean and disinfect.
- LCD display that is easy/ convenient to use and read.
- Display comes with lighting, allows clear visibility in the night.
- Small, handy to use and mobile.
- Suitable for clinical use, on Hygiene Pulp™ products.

### PRODUCT SPECIFICATIONS

REF	Dimensions (L x W x H)	Average Weight/Pc	Colour
05901	200 mm X 165 mm X 19 mm	~388g	3 colour print, White (base), black and orange

### Composition

- ABS plastic (Bottom)
- 3mm tempered glass (Top)
- AAA battery
- LCD display

### Country of Origin

- China

### Cautions

- Do not use any chemical/abrasive cleaning agents to clean the scale, other than the glass surface.
- Handle with care. Do not drop the scale.
- This scale is not waterproof. Do not immerse in water.
- Do remove battery if the scale is not being used for a long period of time.
- Do not frequently overload the scale.

### Care and Cleaning

- To clean: Spray soap solution on the glass surface. Then use another cloth to wipe off.
- Should there be spillage on the scale, wipe it off immediately. The scale is not waterproof.
- To disinfect: Use disinfection wipes or spray on main glass surface. For the rest of the surface use disinfection wipes.
- Do not use any chemical/abrasive cleaning agents to clean the scale, other than the glass surface.

### Shelf Life/ Warranty Shelf Life

- 8 years (Shelf life)
- 2 years warranty

### Selected Standards & Guidelines

- CE, RoHs, EMC, LFGB, PAHS, Safety Data Sheet (AAA Battery)
- Fragmentation test (tempered glass)


### Recommended Storage Condition

- Please keep it dry and at room temperature.
- Avoid direct sunlight.

### Proposed Disposal Method

- Do not treat it as household waste and throw into garbage bag.
- Send to applicable collection point for the recycling of electrical and electronic equipment.

### Product Symbols

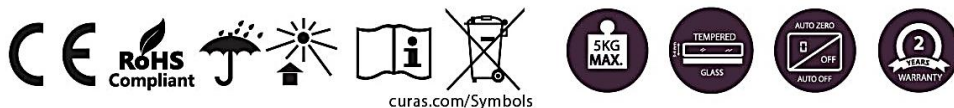

### PACKAGING INFO

REF	EAN-13 Barcode	Carton	Gross Weight
05901	5051400059018	20 pcs / carton	~10.4 kg

Components / Material Properties			Recommendations on Waste Handling			
Packaging	Material	Contains PVC	Incinerate	Landfill	Recycle	Reuse
Inner Box	Thin corrugated	--	√	(√)	√	
Outer Carton	Corrugated Carton	--	(√)	(√)	√	
Transport Pallet	Wood	--			(√)	√

√ = Highly recommended

(√) = Recommended

### ADDITIONAL INFO

- HS Code: 8423101000
- Documentation available upon request. Kindly contact: [info@curas.com](mailto:info@curas.com)