

QUICKIE[®]
Live without limits.

JIVE[®]

The ultimate
driving experience



The **JIVE**[®] family
of high performance
powerchairs



JIVE M



JIVE R



JIVE F



The first time you experience a Quickie Jive, it all becomes clear. It's a powerchair that effortlessly blends indoor agility with outdoor high performance.

Fully customisable and available in a choice of front, rear and mid wheel drive bases, you're sure to find the perfect Jive to meet your needs. Here are just a few of the options available to you...

Choose your seating system

JAY Comfort



JAY clinical seating



Recaro™ and other third-party seating



Choose your powered module

Lift/tilt modules



30° or 50° centre of gravity tilt module



Choose your drive base

Rear wheel drive



Mid wheel drive



Front wheel drive



Enjoy the all-inclusive package

Jive comes fully loaded with some high-end features as standard, including:

- JAY Comfort optimum comfort cushion
- JAY Comfort backrest
- Headrest
- Lap belt
- A choice of angle adjustable footplates or centre mounted legrest
- Parallel swing-away controller
- 4mph (6kph) or 6mph (10kph) high torque (4 pole) motors
- 60Ah batteries
- Solid puncture proof wheels and casters
- LED lights
- Rear-view mirror

JIVE® M

Indoors or out, it's the
best of both worlds



With an ultra-small turning circle, excellent indoor and outdoor handling and superb suspension through SpiderTrac™, the mid-wheel drive option is ideal for those who want a great all-round powerchair.




Ideal for indoors

With Jive M the centre of gravity is located directly beneath you, which means it can turn on its own axis. This ensures an ultra-small turning circle for excellent manoeuvrability around the home.

Offering an intuitive drive, Jive M is ideal for those who may have difficulty driving a powered wheelchair or using switches or head controls.

Navigating through doorways, into lifts and around confined spaces is easy with the Jive M's narrow width of just 62cm (24").

A low seat height of 43cm (17") (with 30° and 50° tilt module options) helps with access under tables and transportation within vehicles, whilst a 90° leg angle ensures good leg positioning.

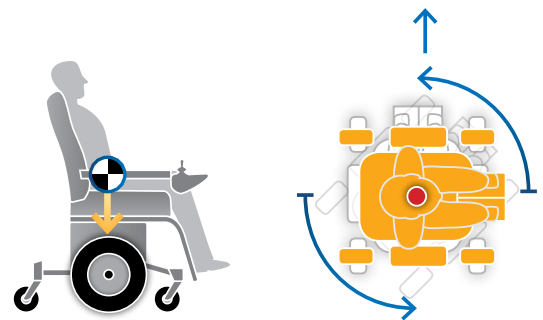



Perfect for outdoors

The SpiderTrac™ gas suspension system absorbs all those bumps as you travel along - delivering possibly the smoothest ride on the market - even on uneven terrain.

SpiderTrac™ delivers gradual, smooth and stable descent down kerbs of up to 10cm (4").

60Ah or 70Ah batteries plus 6mph (10kph) high torque motors mean you'll effortlessly climb slopes of up to 10° and power over rough, rugged terrain with ease.

The seating position also adds more weight to the drive wheel, perfect for increasing traction and stability, and decreasing sliding when outdoors.

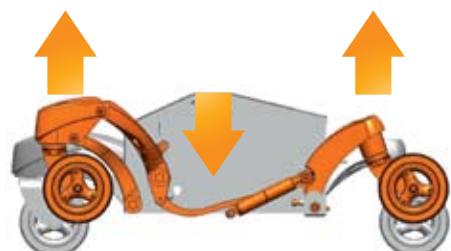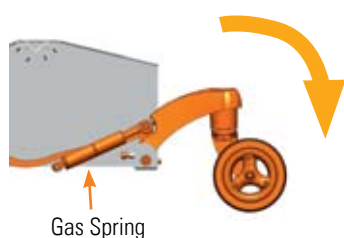




How SpiderTrac™ suspension works:

SpiderTrac™ suspension is a gas suspension system used to dampen the impact without adding pressure to the castors. By maximising the articulation range of the front and rear casters, the Jive M has a spider-like ability to crawl smoothly onto and over steep transitions - so when you're ascending or descending a kerb, the gas dampener cushions the impact of the caster wheel. This ensures not only super-smooth transitions without jolts, jarring and forward pitching but also keeps you seated in a horizontal position.

Because the suspension on the front and rear casters is linked, they move in an opposite direction to compensate for drive wheel movement. So when you're transitioning a slope (or driving from the floor up a car ramp) every wheel maintains contact with the ground, maximising traction. The gas suspension is also used to absorb shock for both casters, providing a more comfortable ride on uneven terrain. You really have to try it to believe it.

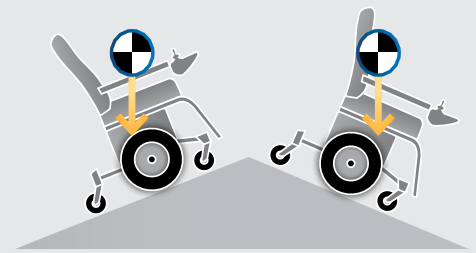

A great all-rounder

- Clean, slim and sleek design
- Controller with optional R-net modules (including Bluetooth) mounted into the base
- Increased foot space between the front caster wheels
- Adjustable centre of gravity on all bases
- Independent adjustable all-wheel suspension
- Speeds of 4mph (6kph), 6mph (10kph) and 8mph (13kph)

Easy to care for

The Jive M is a low-maintenance option which is easy to service. With easy access to batteries and mechanical components, the powerchair also features an automatic snap-in seat stay mechanism.





Powered options

As you would expect, Jive M offers a full range of powered options

- Powered centre of gravity tilt of 0 to 30° (with a -5 to 25° option)
- Centre of gravity tilt of 0 to 50° to optimise pressure release
- Powered lift of 30cm (11.8")
- Powered lift and tilt (30° tilt and 30cm lift)
- Powered reclining back (-10 to 40°)
- Powered anti-shear biomechanical back recline (-10 to 40°)
- Powered legrests
 - Elevating and articulating
 - Centre mount



Jive F and Jive R:
High-end performance
to suit YOU.



Choosing your drive base is a very personal decision. It can be down to situation, or simply just preference. That's why Jive is also available in both front and rear wheel drive options. Both share some very similar features, such as their narrow 64cm (25") width, but there are some unique differences.



Jive R

If you like the familiarity of rear wheel drive, opt for the Jive R – it's intuitive and geared for speed without compromising on stability or control. Indoors the low seat to floor height of 43cm (17") (with tilt) ensures easy under-table access.

Jive F

If you prefer a little extra muscle when climbing those kerbs or tackling off-road terrain, then the Jive F may be for you. With a leaner turning radius of 70cm (27.5") (compared to 100cm (40") on the Jive R), you'll be able to achieve tighter cornering for greater manoeuvrability indoors.

Whatever drive base you choose, you'll enjoy an incredibly responsive driving experience.

Great outdoor performance

- 43cm (17") seat to floor height with 30° and 50° tilt module (RWD) option
- 60 or 70Ah batteries
- 4mph (6kph), 6mph (10kph) and 8mph (13kph) (RWD option only) speed packages
- 90A or optional 120A R-net controller
- Overall width of 64cm (25")

Clean and accessible design

- Clean top shroud design for easy and simple battery/controller access
- New auto-lock seat stay mechanism
- New shroud design that contains R-net modules in the base

Power options

- Powered centre of gravity tilt of 0 to 30° (with a -5 to 25° option)
- Powered lift of 30cm (11.8")
- Powered lift and tilt (30° tilt and 30cm (11.8") lift)
- Powered reclining back (-10 to 40°)
- Powered anti-shear biomechanical back recline (-10 to 40°)
- Power legrests
 - Elevating and articulating
 - Centre mount



Powered centre mount legrest



Seating systems

Sleek, highly adjustable, colour-coded to match your chair and supremely simple. Adapting seating systems and accessories on your Jive couldn't be easier with its versatile seat frame

Designed to fit you, the Jive's seat frame provides simple adjustment of the armrest and seat width/depth (from 40 to 56cm (15.7 to 22"). Not only are adjustments effortless, but the bolts are completely seamless, maintaining the sleek appearance of your seat frame.

Seat frame

- Seat frame offers adjustability via the armrest from 40 to 56cm (15.7 to 22"). Standard cushions are available in 40, 44, 48, 52 and 56cm (15.7, 17.3, 19, 20.5 and 22"). Any additional sizes are available through Built-4-Me.
- Various armrest styles are available, including a multi-adjustable reclining armrest.



Width and depth adjustable armrest

Adding accessories

- Accessories such as oxygen bottle holders or specialty controls are just as easy to mount by using fixation clamps that fit securely onto the seat frame's oval tubing – it's that easy.



Adjustable push handles are just one of the accessories that can be easily fitted into the oval tubing frame.

Oval tubing frame





Jay Comfort seating

This comes as standard and is a range of cushions and back options which provide great support.



Jay Clinical seating option

If a higher level of positioning or skin protection is required, the Quickie Jive seat frame can be quickly adapted to fit a wide range of Jay clinical cushions and backs.



The J3 back will fit into the seat frame and can be easily and quickly removed.

Custom seating

Seat and back attachments are designed for easy adaptation to easily accommodate custom mouldings or third-party seating systems.



Recaro™ seating

Recaro™ seating is available with any of the three drive bases.

Backrest mounting brackets allow simple attachment of custom made seating.



Control systems

With Jive, control is all about simplicity. From the slimline, accessible design of the joystick, to adjusting the seating position or speed, easy control is at your fingertips.

You can further customise your Jive by expanding the control system with R-net. Not only can you control your powered seating, you can also enjoy the increased power and performance this delivers on uneven terrain and when climbing kerbs. Plus you can personalise your Jive to suit you by assigning commonly used commands and settings to a single button.

For conversations on the move, keep tweets and status updates flowing by using the Jive's Bluetooth option to control your smartphone or tablet.

But online social interaction is just the beginning – the environment around you will react to your every command with the Jive's Omni controller. Getting ready to watch a film? Dim the lights, draw the curtains, select play and turn up the volume – all without your hand leaving the joystick. And, if a joystick's not suitable you can even customise the way you control your Jive with a variety of speciality controls.

VR2 non-expandable controls (standard)

If you just need a high performance base with a basic controller function the VR2 standard controller is for you. Dependable, highly flexible and simple to use, it provides 90A motor output for enough torque to manage smooth or uneven terrains.



VR2 std joystick



Dual control option

Optional R-net (120Ah) expandable controls

For more functionality and higher torque choose the modular R-net control system.



LED joystick



LCD joystick



Full colour Omni + control interface

- Features up to six actuator controls via joystick or external programmable buttons. Provides simple programming options, outstanding drive algorithms and proven reliability.
- Bluetooth mouse control through the joystick for Windows PCs. Also valid for Android/Windows based tablets and smartphones through additional software.
- Allows unique programming of profiles for easy transitions between environments.
- Free customised Assignable Button functionality for convenient control of common functions.
- OMNI special controls interface with infrared control option (programming database or code learning).



Special controls/input devices

If limited hand usage is an issue but head movement is good there are various different types of controls and switches that can be mounted in the headrest which can be activated by head movements. In addition to the special controls offered on the Jive order form, customised options are also available through Built-4-Me.



Head control



Switch control



Powered swing-away 0-Way proportional joystick



Powered swing-away chin control

0-Way proportional joystick

Use as a chin or hand control. The potentiometer has an adjustable operating force from 10 to 50g and operates with virtually no deflection.

WebIQ online configurator

With our WebIQ online configurator and download platform, Quickie powerchairs can be re-configured to meet new requirements, reprogrammed with default settings and user languages can be changed. Personalised settings can also be uploaded to the chair's serial number identity.



Transit approved

The Jive family of powerchairs have been successfully crash tested and fulfil the performance requirements for ISO 7176-19.

They are therefore suitable to be used as a seat in a vehicle for the transportation of wheelchair users.

Customising

The Jive can be further customised through Built-4-Me. Whether the request is for a clinical need or simply to give the powerchair a look which makes it stand out from the crowd, Built-4-Me can handle it. Simply ask your Quickie provider to contact the Built-4-Me team with your specific request, and our experts will then assess and discuss any risk factors involved before designing and manufacturing a product to meet your exact needs.

Examples of customising

These are just a few examples of Built-4-Me customised options for Jive.



Customised colours for bases, also available for Jive F and R frames.

Seating in different materials and colours.



For more information on the custom built option please visit www.customwheelchairs.eu or contact your Sunrise representative.

Control adaptations

These are just a few examples of possible customisations for Jive controls.



Chin control powered swing-away with individual position of switches



Tray mounted flip-up joystick



Centre mount joystick with actuator box



Switches and joystick mounted in tray

Sport adaptations such as the football or hockey kit

Choose your team colours, select the options and controls you need, add the protection and drive profiles to aid your performance and you could have a Jive Sport built specifically for your recreational needs.

Specialist sports options include:





- Rear ram protection
- Front ram protection (including footrest bar)
- Anti-pitch blocking

And for the ultimate in flexibility the anti-pitch blocking and bumpers can be easily removed to convert the Jive into an everyday powerchair in just five minutes.




Technical Specifications:


F&R: 40 - 56cm [15.5 - 22"]
M: 40 - 56cm [15.5 - 22"]

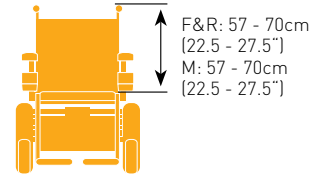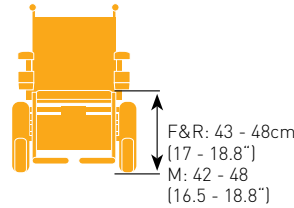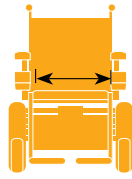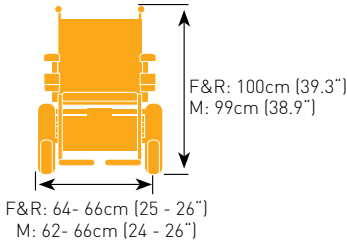

F&R: 40 - 56cm [15.5 - 22"]
M: 40 - 56cm [15.5 - 22"]

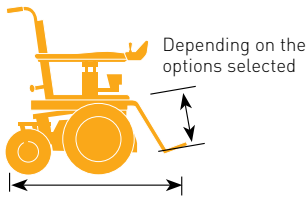

F&R: 118.5 - 121.5cm [46.5 - 47.5"]
M: 107 - 113cm [42 - 44.5"]

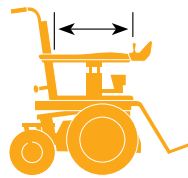

F&R: 18 - 28cm [7 - 11"]
M: 18 - 28cm [7 - 11"]

F&R: -10° - +40°
M: -10° - +40°



F&R: 0° - +30°
M: 0° - +30° / 0° - +50°

F&R: 88cm [34.5"]
M: 95cm [37.5"]

Speed:	Jive R&M: 6, 10, 13kph	Jive F 6, 10kph	Total weight (with batteries):	145 - 187kg (319 - 412lbs)
Batteries:	60Ah/ 72Ah		Max kerb height:	10cm (4")
Max. range:	Approx. 25 miles (40km)		Max safe slope:	F: 6° (11%), R & M: 10° (18%)
Max user weight:	25 stone (160kg)	Minimum turning radius:	M: 56cm (22"), R: 100cm (39.3"), F: 70cm (27.5")	

Find out more today! Visit www.QuickiePower.com/Jive for more information, to find your local dealer, or to request a demonstration.



Get SOCIAL with QUICKIE

Keep up-to-date with the latest news on the Quickie Jive and other Quickie products through our Facebook page:

www.facebook.com/QuickieWheelchairsUK



Sunrise Medical Limited
Thorns Road
Brierley Hill
West Midlands
DY5 2BQ
Tel: +44 (0) 845 605 66 88
Fax: +44 (0) 845 605 66 89
www.SunriseMedical.com